CanalPro™ Apex Locator Compact

3	Gebrauchsanweisung
17	Instructions for use
31	Mode d'emploi
45	Istruzioni per l'uso
59	Instrucciones para el uso
73	Instruções de utilização
87	Gebruiksaanwijzing
101	Brugsanvisning
115	Bruksanvisning
129	Bruksanvisning
143	Käyttöohje
157	Οδηγίες χρήσης
171	Instrucțiuni de utilizare
185	Instrukcja użycia
199	Navodila za uporabo
213	Naudojimo instrukcijos
227	Kasutusjuhised
241	Návod k použití
255	Használati utasítás
269	Návod na použitie
283	Инструкции за употреба
297	Lietošanas instrukcija
311	Kullanma talimatları
325	Upute za upotrebu

CanalPro™ Apex Locator Compact

Inhalt

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 6.1. Packaging box content..... 20

 6.2. Getting started 20

 6.3. Installing / replacing the batteries 21

 6.4. Device operation..... 22

 6.5. Cable connection test..... 22

 6.6. Apex localization..... 23

 6.7. Audio feedback..... 25

 6.8. Demo mode..... 25

 6.9. Automatic shutdown 25

7. Maintenance, cleaning and sterilization..... 25

8. Troubleshooting 27

9. Warranty 27

10. Disclaimer..... 28

11. Certification 28

12. Technical specifications 28

13. Standard symbols..... 29

Directions for use

■ Introduction

Congratulations, you have made a good choice with CanalPro Apex Locator Compact. The ergonomic user-friendly design of CanalPro Apex Locator Compact makes it a preferred choice in modern dentistry.

Utilizing the advantages of well proven patented technology, CanalPro Apex Locator Compact is a precise, user-friendly, affordable apex locator suitable both for skilled endodontic experts and for general practitioners performing root canal treatments. The custom color graphic display of the CanalPro Apex Locator Compact helps to achieve optimal performance during determination of the working length.

1. Indications for use

The CanalPro Apex Locator Compact is an electronic device used for apex localization and working length determination during root canal treatment. The device enables you to obtain correct results in canals under various conditions – dry and wet.

2. Contraindications

The CanalPro Apex Locator Compact is not recommended for use in patients with a pacemaker or other implanted electrical devices.

3. Warnings

This product must only be used in hospital environments, clinics or dental offices by qualified dental personnel.

4. Precautions

- Do not use the CanalPro Apex Locator Compact near devices emitting electromagnetic noise such as fluorescent lamps, film viewers, ultrasonic devices, etc. Cellular phones, remote controls or other devices generating electromagnetic waves may interfere with the operation of the CanalPro Apex Locator Compact. Such devices should be switched off.
- During operation, protect the CanalPro Apex Locator Compact from possible spillage of liquids.
- Do not use the CanalPro Apex Locator Compact in the presence of flammable materials.
- The CanalPro Apex Locator Compact should only be used with the manufacturer's original accessories.
- In order to prevent the transmission of infections it is highly recommended to use a rubber dam system during the endodontic procedure.
- To ensure that short circuits do not impair the measurements, avoid any contact between metal crowns, bridges or amalgam fillings and the file or the lip clip.
- Make sure that the canal is wet enough to provide electrical conductivity.
- Ensure that the file does not touch other instruments.
- Avoid excessive liquids inside the tooth cavity to prevent overflow and incorrect measurements. Teeth with open apices may give imprecise results.
- Usage of NaOCl with concentrations higher than 5 % may result in reduced accuracy of the measurements.
- The CanalPro Apex Locator Compact may not be able to provide correct measurements in all conditions. In all cases, it is recommended to take an X-ray prior to using the device and to compare the results obtained by both methods.
- For your own safety, please use personal protection wear (gloves, mask).
- The CanalPro Apex Locator Compact should be stored at normal temperature (< 60°C) and humidity conditions.

5. Adverse reactions

None.

6. Step-by-step instructions

■ 6.1. Packaging box content

Check the content before use:

- CanalPro Apex Locator Compact apex locator 1 pc.
- Base for device 1 pc.
- Lithium coin battery CR2032 2 pc.
- Measuring cable 1 pc.
- Lip clip 2 pcs.
- File clip 2 pcs.
- Touch probe 1 pc.
- Screw driver 1 pc.
- User manual 1 pc.

■ 6.2. Getting started

3 pushbuttons are located on the front panel:

Sound

Mode

On/Off

disconnect the charger from the device, if connected

The connection port for the measuring cable is located on the right side of the device. The battery compartment is located on the rear side of the device. The device can be used with or without the base.

Connection port for measuring cable

■ 6.3. Installing / replacing the batteries

The CanalPro Apex Locator Compact is powered by 2 standard Lithium coin cell batteries, type CR2032.

- ⚠ Prior to battery replacement, the device must be switched off.
- ⚠ Before replacing the batteries, the measuring cable should be disconnected from the patient and removed from the CanalPro Apex Locator Compact device.

6.3.1. To install / replace the batteries, release the screw and remove the battery compartment cover on the rear side of the device (Fig. 1a).

- Remove old batteries from battery holders using a small screwdriver as shown in Fig. 1b.
- Insert each of the new batteries into the battery holder, slide forward and down until fully seated in the holder.
- Close the battery compartment and secure the cover with the screw.

- ⚠ Ensure correct polarity (the "+" sign on the battery must face upwards).

■ Figure 1a

Figure 1b

6.3.2. When the batteries are depleted, a flashing Low Battery indicator will appear in the upper right corner of the screen (Fig. 2). The CanalPro Apex Locator Compact will continue normal operation even with low batteries, but the batteries should be replaced before becoming flat.

6.3.3. If the batteries are flat, the flashing Low Battery indicator will be accompanied by 3 beeps and the device switches off automatically.

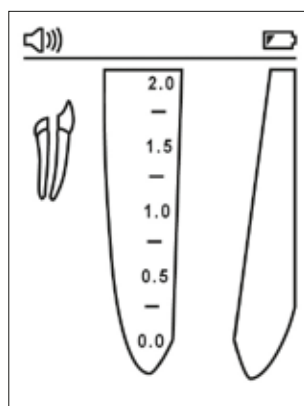
Figure 2

■ 6.4. Device operation

Prior to first use, it is recommended to sterilize the Lip Clip, the Touch Probe and the File Clip. Please refer to section 7 for further information regarding CanalPro Apex Locator Compact maintenance.

- 6.4.1. Switch the device on by pressing the On/Off button . The main screen is displayed. The status bar with cable, battery and sound level indicators is located in the upper part of the display.
- 6.4.2. Insert the measuring cable into the connection port of the device and make sure that the cable icon appears on the status bar (Fig. 3).
- 6.4.3. Connect the lip clip and the File Clip/touch probe to the measuring cable.
- 6.4.4. Attach the Lip Clip to the patient.
- 6.4.5. Insert endodontic file into the canal and connect the file clip to the file (to ensure precise measurements, the file size should be adjusted to the canal diameter). Alternatively, use the touch probe to provide electrical contact with the endodontic file.

The file image inside the small tooth icon will appear accompanied by a double beep audio signal. The absence of the file image and of audio feedback indicate a faulty connection. Check your connections, clean your file clip/touch probe contact, moisten the canal, if necessary, and start again.

No other adjustments are required before starting measurements.

■ 6.5. Cable connection test

A connection test feature is included in the CanalPro Apex Locator Compact in order to check the cables:

- Insert the measuring cable into the connection port and connect the file clip and the lip clip to the cable.
- Verify that the file clip contact and the lip clip are clean and establish a connection.
- The “Connection” symbol  should appear on the status bar of the display, indicating proper connection (Fig. 4).
- If the symbol does not appear, the measuring cable or the file clip should be replaced.

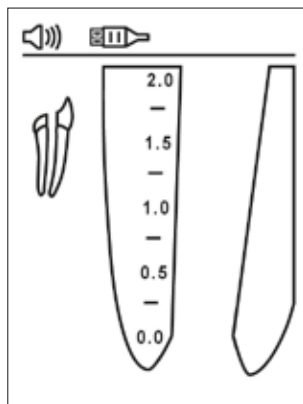
Figure 3

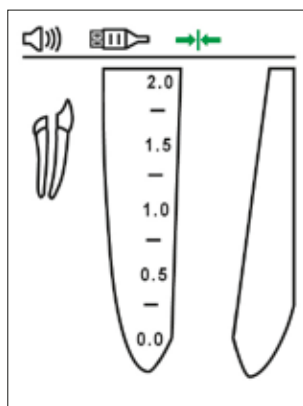
Figure 4

■ 6.6. Apex localization

Slowly introduce the file into the canal. Movement of the file inside the canal is reflected on the screen by green arrows Down and Up (Figs. 5 and 6).

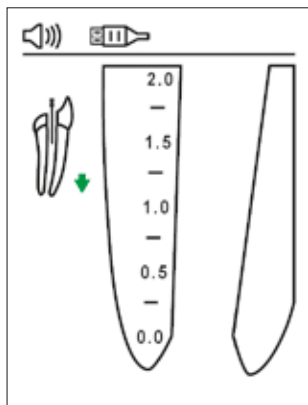
Figure 5

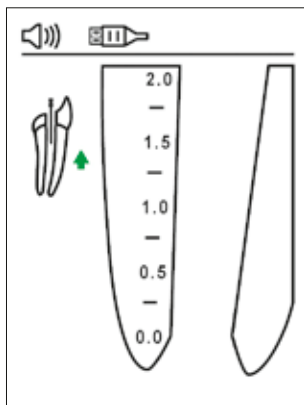
Figure 6

Along the pre-apical zone, a zoomed view of file progression in the canal is shown on the enlarged root image in the right side of the display, and the numerical value appears on the left side of the display under the tooth image.

The CanalPro Apex Locator Compact provides audio feedback of file progression as a series of beeps with variable rate.

Beginning of pre-apical zone

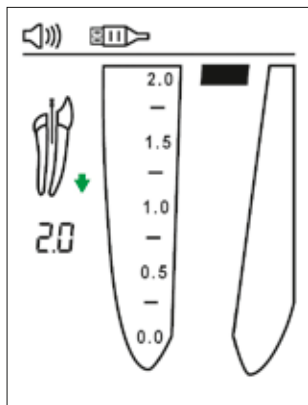
Figure 7

Pre-apical zone

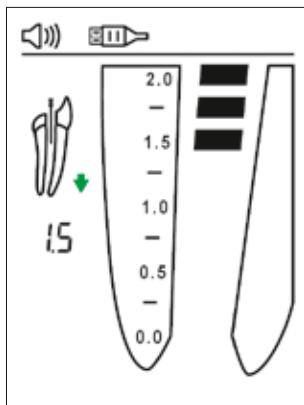
Figure 8

6.6.1. Apical zone

File position in the apical zone is visually indicated by color segments inside the root canal image and numerical readings 1.0, 0.75, 0.5, 0.25 and 0.0 under the tooth image in the left part of the display.

File progression in the apical zone is accompanied by a series of beeps with variable rate. When the apical position is reached (red bar at the mark 0.0), a constant tone is emitted.

Beginning of the apical zone

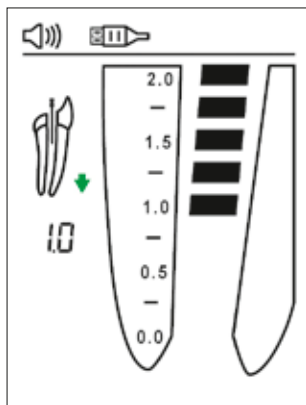
Figure 9

Mid-apical zone

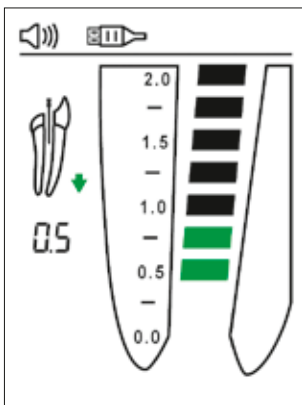
Figure 10

Apical position

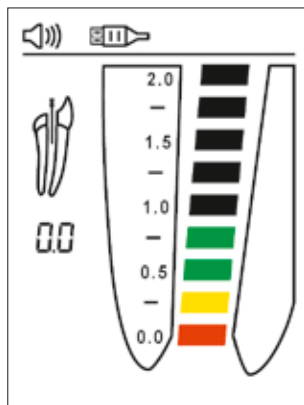
Figure 11

The apical position indicated by the CanalPro Apex Locator Compact (red bar at the mark 0.0) corresponds to the file tip position at the minor apical foramen.

Determining the working length is a matter of the dentist's professional judgment. In most cases, subtraction of 0.5 mm from the apical length provides a clinically acceptable working length. Nevertheless, in each case the dentist should define a proper working length based on his experience, apex locator readings, radiographs and other available data.

6.6.2. Over-instrumentation

A red "Blood drop" and warning sound indicate that the file has passed the apex. An "O.A" reading appears in the left part of the screen under the tooth image (Fig. 12).

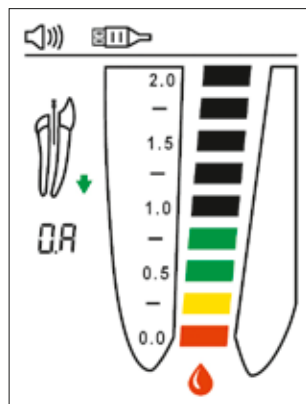
Figure 12

■ 6.7. Audio feedback

The CanalPro Apex Locator Compact is equipped with an audio indicator which enables monitoring the progression of the file within the canal in the pre-apical and apical zones.

The sound volume can be adjusted to one of four levels: mute, low, normal and high, by successively pressing the volume button .

■ 6.8. Demo mode

The built-in Demo mode is available to demonstrate operation of the device and to improve the learning curve of the operator.

- 6.8.1. Disconnect the measuring cable from the CanalPro Apex Locator Compact, if connected, and switch the device on.
- 6.8.2. To start the Demo mode, press and hold the Mode button  for about 1 sec. until the second beep sounds and "Demo" indication appears on the screen.
- 6.8.3. During the Demo cycle the operating sequence of the device is shown on the screen. Press Mode button  to pause the simulation; press  button again to resume.
- 6.8.4. When the Demo cycle is completed, it is repeated automatically until interrupted by the operator.
- 6.8.5. To exit the Demo mode, press and hold the Mode button  for about 1 sec. until a beep sounds.

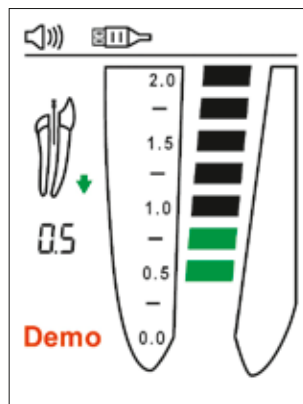
Figure 13

Note:

If the measuring cable is plugged into the connection port of the device during the Demo cycle, the CanalPro Apex Locator Compact will exit the Demo mode and automatically switch to regular operation mode.

■ 6.9. Automatic shutdown

The CanalPro Apex Locator Compact automatically shuts down after 2 minutes if not in use. In order to prolong the battery life, it is recommended to switch off the device after completing the measurements by pressing the On/Off button .

7. Maintenance, cleaning and sterilization

■ 7.1. General

- The device does not contain parts which can be serviced by the user. Service and repair should be provided by factory trained service personnel only.
- All objects that were in contact with potentially infectious agents should be cleaned after each use:

Lip Clip, File Clip and Touch Probe should be disinfected and sterilized by autoclaving between treatments. Please follow "Disinfection and sterilization procedure" described in section 7.2.

Measuring cable, the device and the station should be cleaned using a tissue or soft cloth impregnated with aldehyde-free disinfecting and detergent solution (bactericidal and fungicidal).

 The measuring cable cannot be autoclaved.

 Use of agents other than specified above may cause damage to the equipment and its accessories

■ 7.2. Disinfection and sterilization procedure

	Operation	Instructions	Details and Warnings
1	Preparation at the point of use prior to processing	No particular requirements.	
2	Preparation for decontamination/ preparation before cleaning	No particular requirements.	
3	Cleaning: Automated	The accessories are not intended for automated cleaning.	
4	Cleaning: Manual	• Clean the accessories with an adequate brush or towel soaked in a disinfectant solution.	• The file clip should be activated during the cleaning process (pressed and released several times). • After cleaning no visible impurities should remain on the accessories.
5	Disinfection	• Soak the required accessories in a disinfectant solution combined with proteolytic enzyme if possible. • Thoroughly rinse the accessories in flowing water.	• Follow instructions given by the disinfectant manufacturer (concentration, immersion time, etc.). • Do not use disinfectant solution containing aldehyde, phenol or any products which may damage the items.
6	Drying	No particular requirements.	
7	Maintenance, inspection and testing of the accessories	No particular requirements.	
8	Packaging	Pack the devices in sterilization pouches.	• Check the validity period of the pouch given by the manufacturer to determine the shelf life. • Use packaging which is resistant up to a temperature of 141 °C (286 °F).
9	Sterilization	• Steam sterilization at 135 °C (275 °F) for 10 minutes in gravity type autoclave. (Table Top, N type) • Drying time after sterilization – 30 minutes.	• Follow maintenance and operation procedures of the autoclave provided by the manufacturer. • The only sterilization parameters to be used are those that have been validated and provided to the user in this User Manual.
10	Storage	• Keep devices in sterilization packaging in a dry and clean environment.	• Sterility cannot be guaranteed if packaging is open or damaged (check the packaging before using the instruments).

8. Troubleshooting

Please review the suggested solutions before calling customer service.

	Problem	Possible Cause	Solution
1	The device does not switch on by pressing On/Off button.	1. The button may be functioning incorrectly. 2. The battery may be empty. 3. Electronic malfunction.	1. Press the On/OFF button several times. 2. Change the battery. 3. Contact your customer service.
2	The device shuts off while measuring.	1. The device was not used for more than 2 min. 2. The batteries are flat.	1. Turn on the device. 2. Replace the batteries.
3	No sound while measuring.	"Sound control" set at "Mute" level.	Adjust the sound level by pressing the "volume" button  .
4	The device does not show file advance inside the canal	1. Poor electrical contact. 2. Electronic malfunction.	1. Clean contact surfaces of the file, the file clip and the lip clip and perform cable connection test as described in the § 6.5. 2. Contact your customer service.
5	Incorrect measurements	1. Dry/calcified canal. 2. The conductive path for restoration treatment is blocked. 3. Blood or irrigation liquid overflow provides conductive path outside the canal. 4. Deep caries provides conductive path outside the canal. 5. Metal crown or metal filling provide conductive path outside the canal. 6. Perforation. 7. Large lateral canal.	1. Irrigate the canal. 2. Remove old canal filling to open the path, irrigate the canal. 3. Dry the excessive liquid. 4. Block the external conductive path. 5. Avoid contact between the file and metal crown/filling. 6. Remove the file, close the perforation and repeat the measurements, carefully inserting the file into the canal. 7. Try to continue the measurements (gently advance the file towards the apex until regular readings appear).

9. Warranty

The CanalPro Apex Locator Compact is warranted for 24 months from the date of purchase. The accessories (cables, batteries, etc.) are warranted for 6 months from the date of purchase.

The warranty is valid for normal usage conditions. Any damage caused by accident, abuse, misuse, or as a result of service or modification other than by a person authorized by the manufacturer will render the warranty void. The warranty is in lieu of any other warranty expressed or implied.

10. Disclaimer

The manufacturer, its representatives and its dealers shall have no liability or responsibility to customers or any other person or entity with respect to any liability, loss or damage caused or alleged to be caused directly or indirectly by equipment sold or furnished by us, including, but not limited to, any interruption of service, loss of business or anticipatory profits, or consequential damages resulting from the use or operation of the equipment.

The manufacturer reserves the right to implement changes and modifications of the product at any time, to revise this publication and to make changes in the contents hereof without obligation to notify any person of such changes, modifications or revisions.

11. Certification

The CanalPro Apex Locator Compact complies with the following standards: IEC 60601-1 (Safety) and IEC 60601-1-2 (Electromagnetic compatibility), including conducted and radiated immunity tests as specified for equipment of Group 1 Class B.

The CanalPro Apex Locator Compact is covered by the “CE Marking of Conformity” certificate. The device bears the following CE identification mark:

12. Technical specifications

The CanalPro Apex Locator Compact electronic apex locator belongs to the following category of medical devices:

- Internally powered equipment
- Type BF applied parts
- Not suitable for use in the presence of flammable anesthetic mixtures with air, oxygen or nitrous oxide
- Continuous operation
- Ingress of liquids – not protected
- The device is intended for indoor use only
- Environmental conditions during transportation: temperature: -20°C to +60°C (0 to 140°F); relative humidity: 10% to 90%, non-condensing

The CanalPro Apex Locator Compact is intended for use in an electromagnetic environment specified for equipment of Group 1 Class B.

Specifications:

Dimensions:	93 × 54 × 13 mm
Weight:	60 gr.
Type of screen:	Custom Color Graphic LCD
Screen dimensions:	51 × 38 mm
Power source:	2 × CR2032 Lithium coin cell batteries

13. Standard symbols

The following standard symbols are used on the device:

Class II equipment

Type BF applied part

Direct current

Manufacturer

Consult instructions for use

Recycling : PLEASE DO NOT THROW AWAY! This product and all its components must be recycled through your supplier

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Coltene/Whaledent GmbH + Co. KG