YOUR DAILY ROUTINE

Enhancing your workflows through innovations

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This document is ACTEON Manufacturing property.

It is not contractual and is intended to ACTEON® representatives and distribution system.

The data contained in this catalogue are valid at the printing date.

All author rights reserved for all countries.

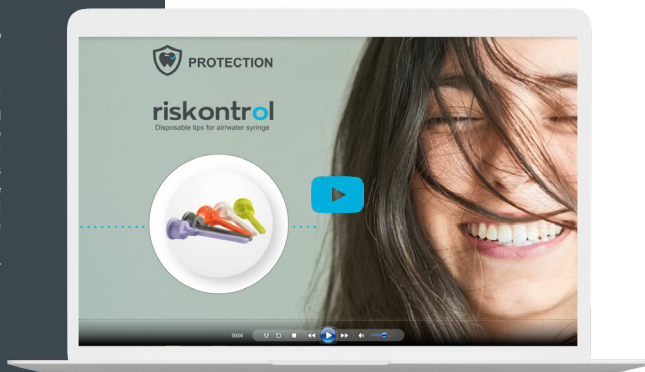
Any entire or part copy made by any way without the writing and previous agreement of ACTEON Manufacturing is unlawful and constitutes a fraudulently copying given rise to penalties.

Adapters: Class I medical devices - CE. RISKONTROL® Classic, RISKONTROL® Art, RISKONTROL® Total Protec: Class I medical devices - CE. RISKONTROL® Sterile: Class Is medical device - CE0459 GMED. Read the instructions in the notice or on the labelling carefully. All information that is essential to ensure correct use of these devices is included in the summary of the product characteristics available on the laboratory's website www.acteongroup.com. For dental professional use only. Not reimbursed by health insurance organisations.

Manufacturer: ACTEON Manufacturing - France. Updated: April 2025.

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*General rules for using
the adapters*

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*Adapters: assembly rules
and drawings*

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*Chairs and syringes brands
and corresponding
adapters*

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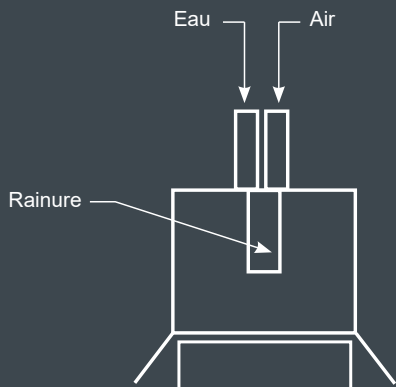
***General rules
for using
the adapters***

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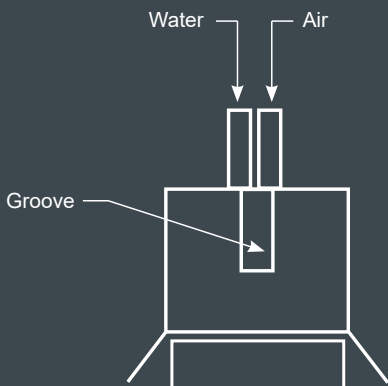
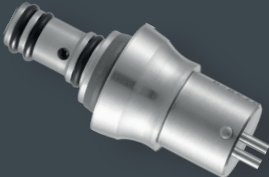
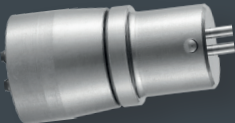
To read imperatively before using the RISKONTROL® adapters

1. Clean/sterilize the adapter before its first use.
2. The adapter must imperatively be set on the syringe according to the assembly drawing before using the RISKONTROL® disposable tips.
ACTEON MANUFACTURING DECLINES ANY RESPONSABILITY IN CASE OF DAMAGES OF:
 - THE SYRINGE DURING THE ASSEMBLING OF THE ADAPTER.
 - THE SYRINGE AND/OR THE ADAPTER IF TIPS OTHER THAN RISKONTROL® TIPS ARE USED.
3. For a good use of the RISKONTROL® adapters, the air and water **pressure** must be **3.5 bars minimum** (about 50 PSI).
4. **IMPORTANT:** Once the adapter in position, check its correct mounting to the syringe before use then verify that the water comes out of the tube the most in the centre of the round part of the adapter (see drawing below).
5. The adapters and syringes with rubber O-rings must be lubricated with silicon oil before the assembling.
6. Between each patient, clean the external part of the adapter with a disinfecting wipe.
Every month, you must remove the adapter from the syringe and lubricate the O-rings.
7. The adapter can be cleaned in an ultrasonic bath nearly every six months after removing it from the syringe and removing the O-rings. For aluminium adapters, use a neutral detergent.
8. The adapter can be sterilized in an autoclave (after removing it from the syringe and removing the O-rings) (132°C - 2.9 bars - 4 minutes or 134°C - 2.2 bars - 18 minutes).
9. The adapter can only be used with the RISKONTROL® disposable tips (sterile or not).
10. The tubes of the adapters can be cleaned out with a scraper Ø 40.
11. Do not discard in household waste. Recycle or dispose of according to current legislations, rather by a waste collector or an approved company.

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Only use these adapters with Riskontrol® disposable tips

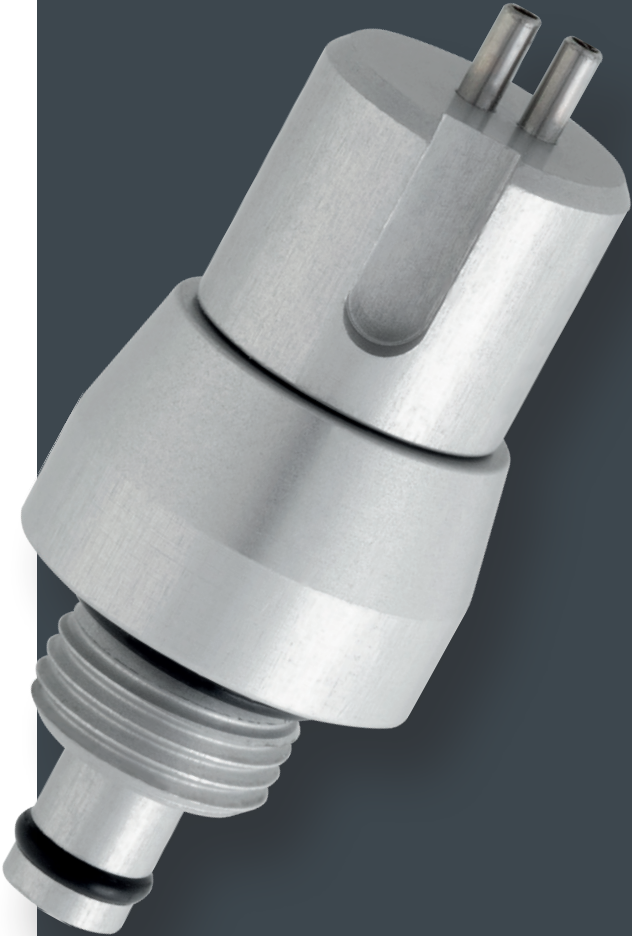
Adaptateur Adapter	Réf. Ref.	Page
Adec Inox	200195	14
Adec New	200205	14
Belmont New	200158	15
Cas Inox	200198	16
Castellini Logo New	200184	17
Castellini V	200163	18
Cefla Inox	200197	19
DCI	200152	20
Dentalez	200164	21
Engle	200194	22
FA C	200167	23
FA C2	200168	24
FA E	200165	25
FA E2	200166	26
FA S1	200182	27
Flex	200169	28
Gnatus Inox	200149	29
Kavo C	200171	30
Kavo DB	200157	31
Kavo E	200172	32
Kavo GB	200156	33
Kavo HB	200153	34
Kavo IB	200 192	35

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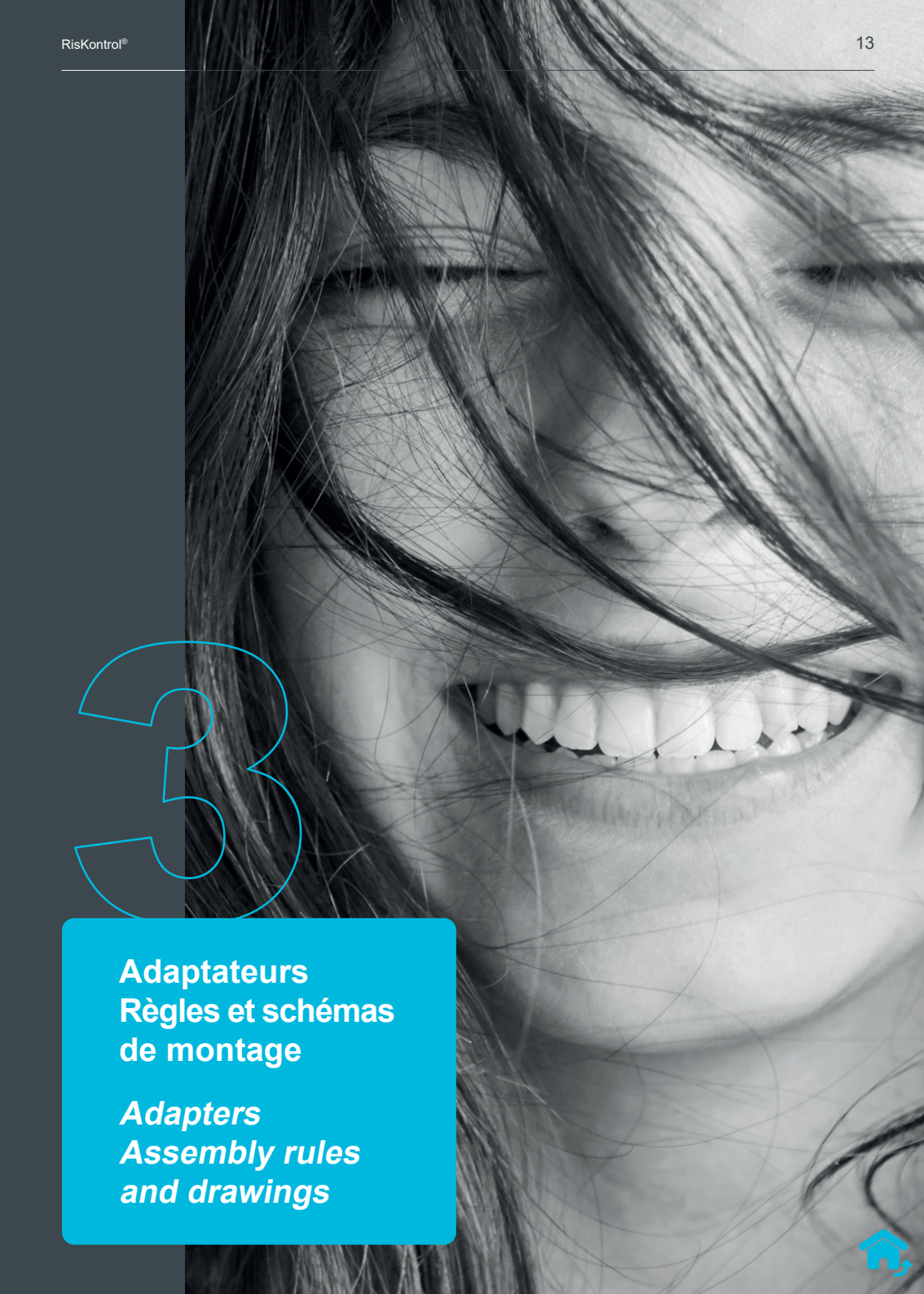
***Adapters
Assembly rules
and drawings***

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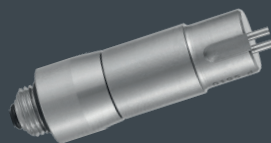
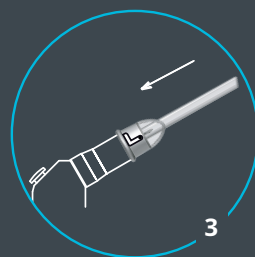
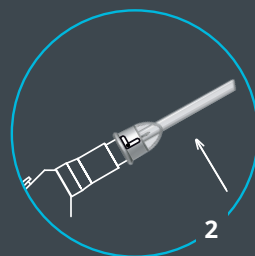
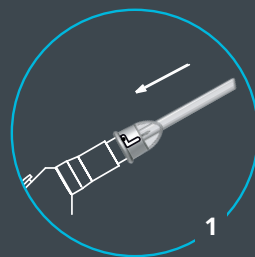
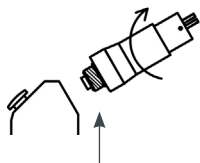
*Pull on the original tip of the
syringe.
Unscrew the nut with an awl key.*

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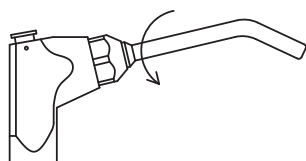
*Check that the small O-ring is at
the end of the adapter.
Lubricate the O-ring and screw
the adapter.*

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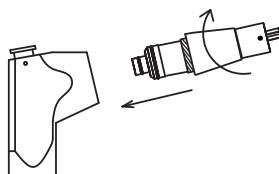
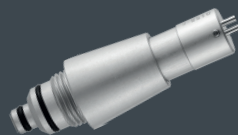
Unscrew the original tip of the syringe.

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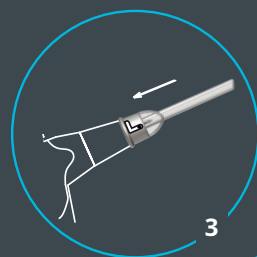
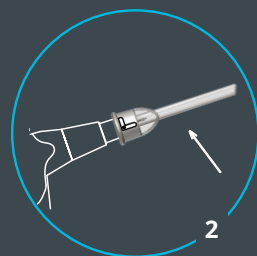
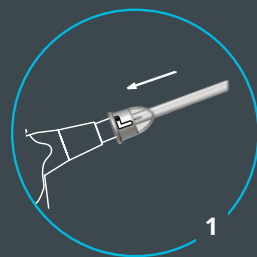
Lubricate the 2 O-rings and push the adapter into the syringe while turning it, then screw the nut from the adapter to the syringe.

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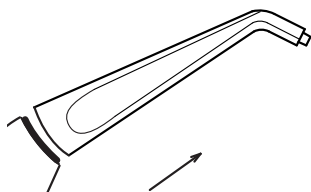
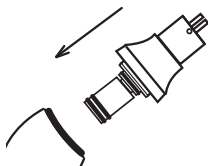
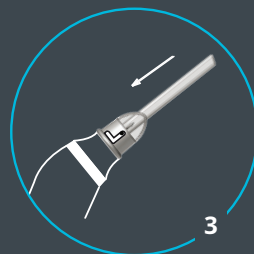
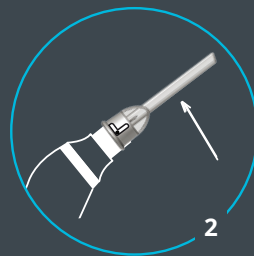
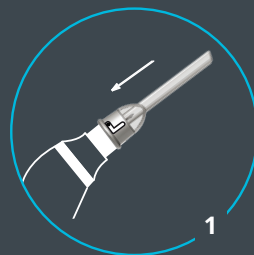
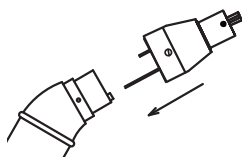
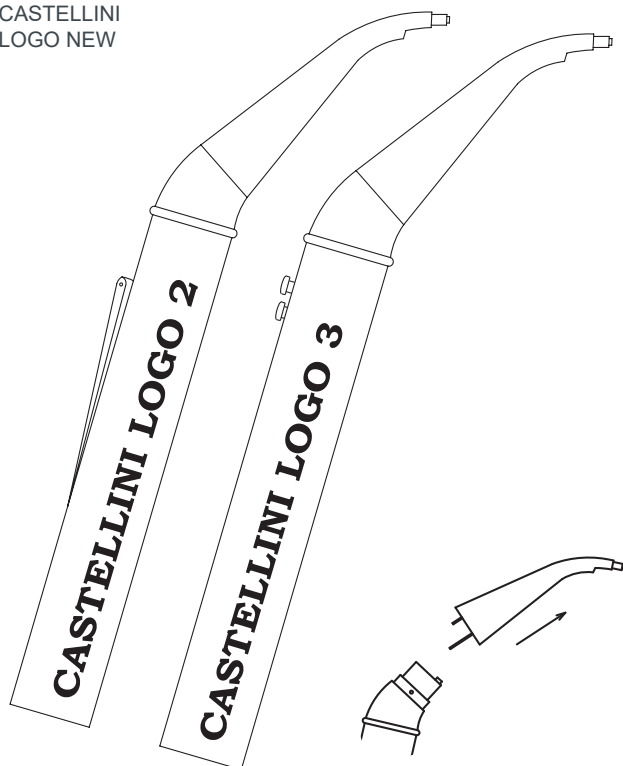
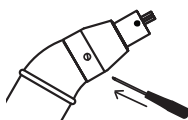
Set up the adapter onto the syringe, putting the 2 tubes into the 2 holes, as for the original tip.

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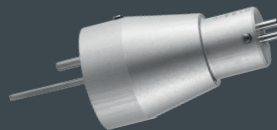
Remove softly the original tip of the syringe by pulling it.

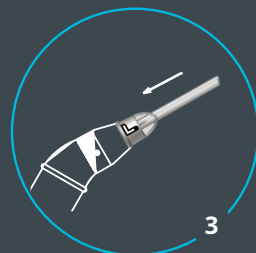
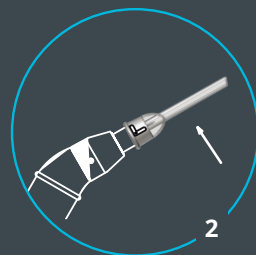
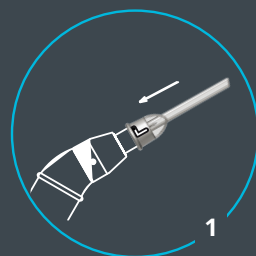
3. Visser la vis de l'adaptateur.
Screw the screw located on the adapter.

Adaptateur sans joint
Adapter without O-ring

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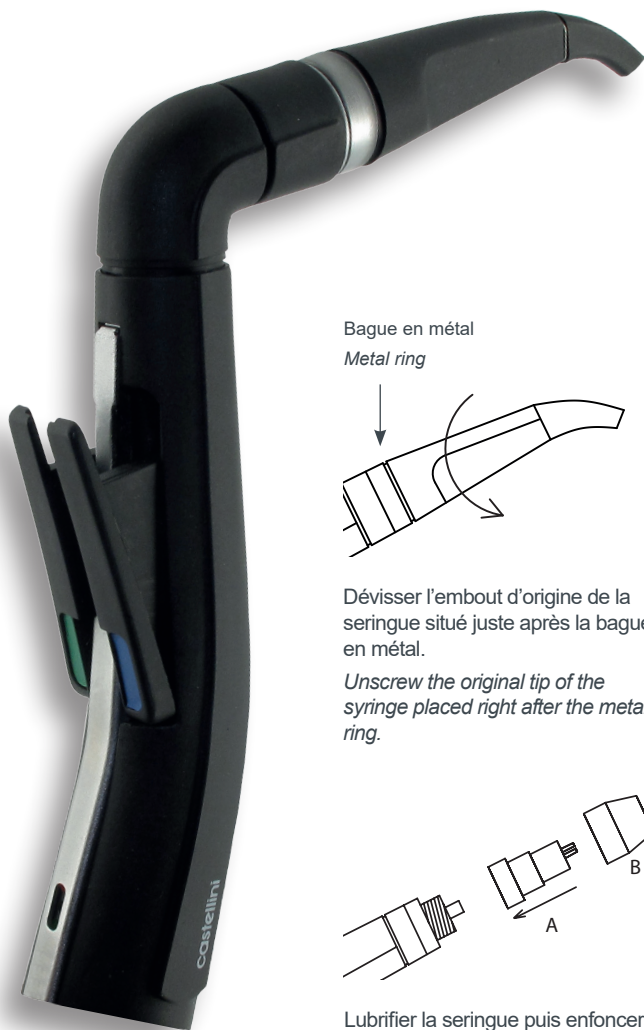
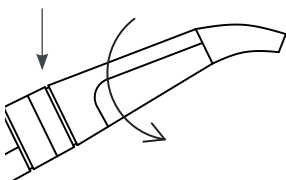
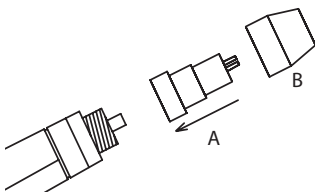
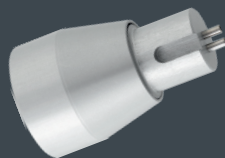
Unscrew the original tip of the syringe placed right after the metal ring.

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Lubricate the syringe and push in part A of the adapter while turning, then screw the part B on the syringe.

Adaptateur sans joint
Adapter without O-ring

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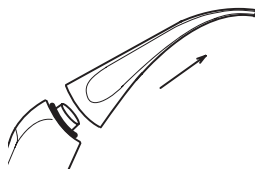
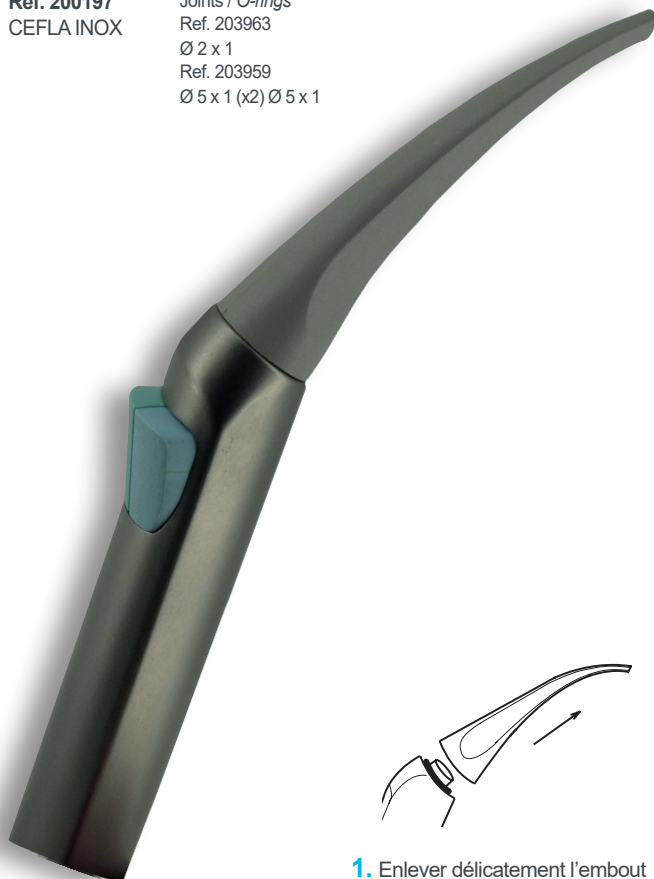
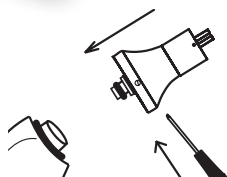
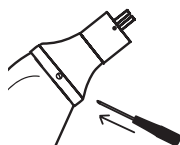
Lubricate the O-rings of the adapter, unscrew the small screw, clip the adapter onto the syringe while turning it.

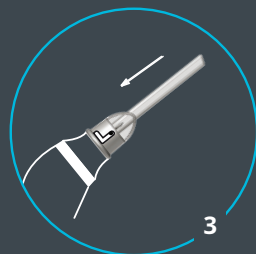
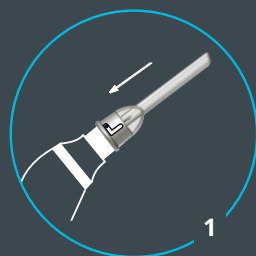
3. Visser la vis de l'adaptateur.
Screw the screw located on the adapter.

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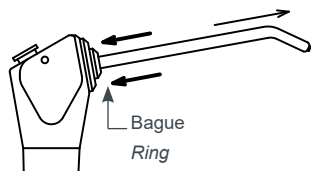
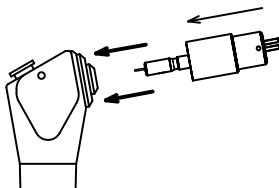
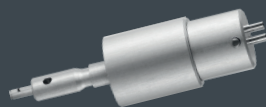
Unclip the original tip of the syringe pressing on the ring around this tip.

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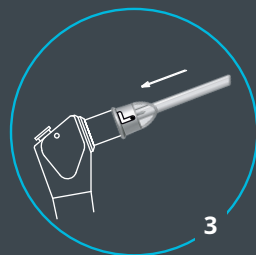
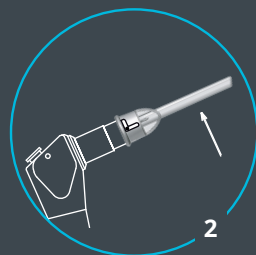
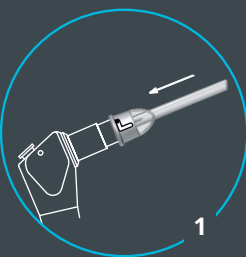
Keep the ring on the syringe before placing the adapter. Clip the adapter pressing on the ring.

Adaptateur sans joint
Adapter without O-ring

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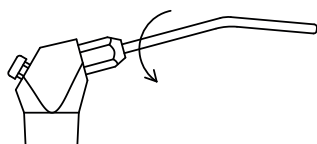
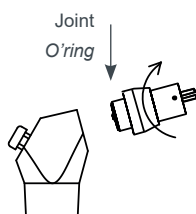
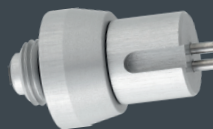
Check that the small O-ring is at the end of the adapter.

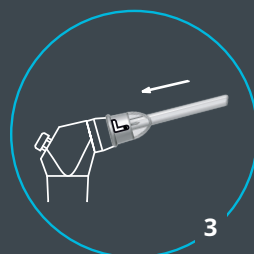
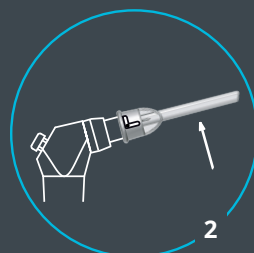
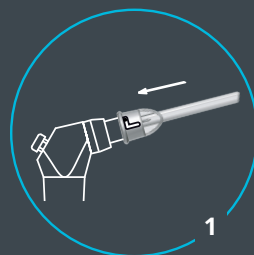
Lubricate the 2 O-rings then screw the adapter.

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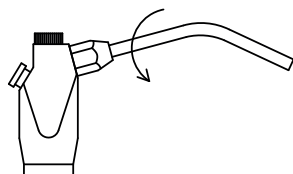
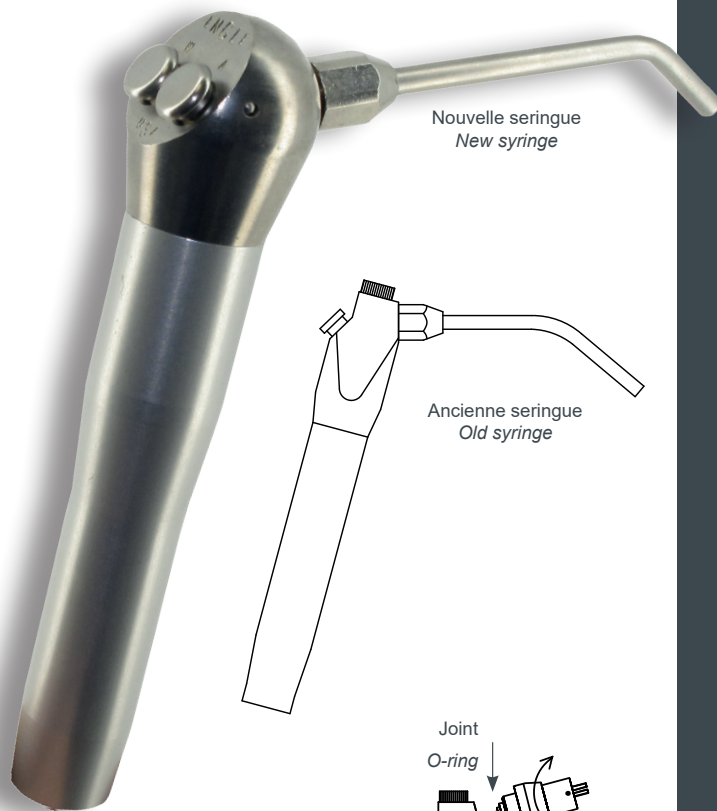
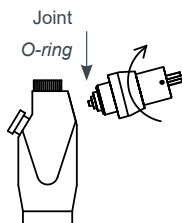
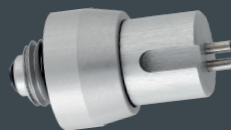
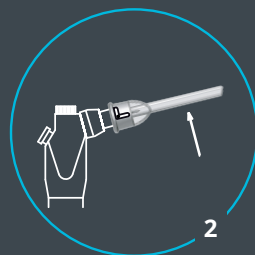
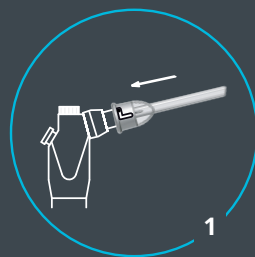
*Check that the small O-ring is at the end of the adapter.
Lubricate the 2 O-rings then screw the adapter.*

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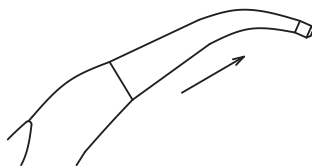
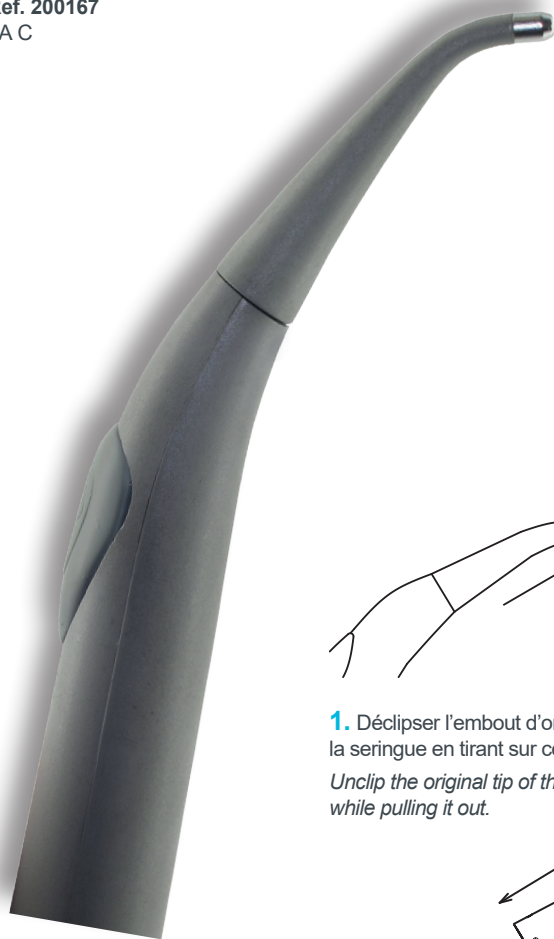
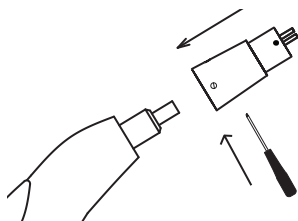
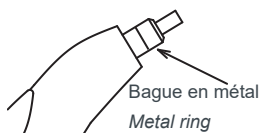
Unclick the original tip of the syringe while pulling it out.

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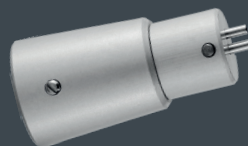
Lubricate the syringe then unscrew the small screw and clip the adapter on the syringe while turning it and screw again the small screw.

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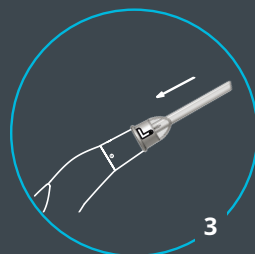
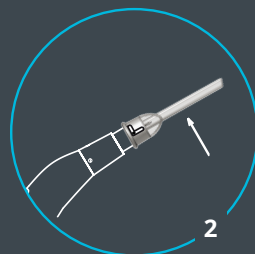
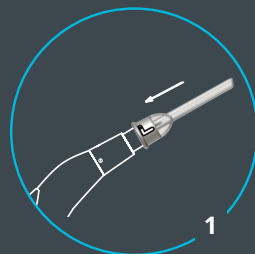
Remove the metal ring (the clip) of the syringe.

Adaptateur sans joint
Adapter without O-ring

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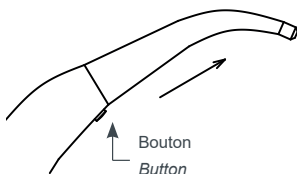
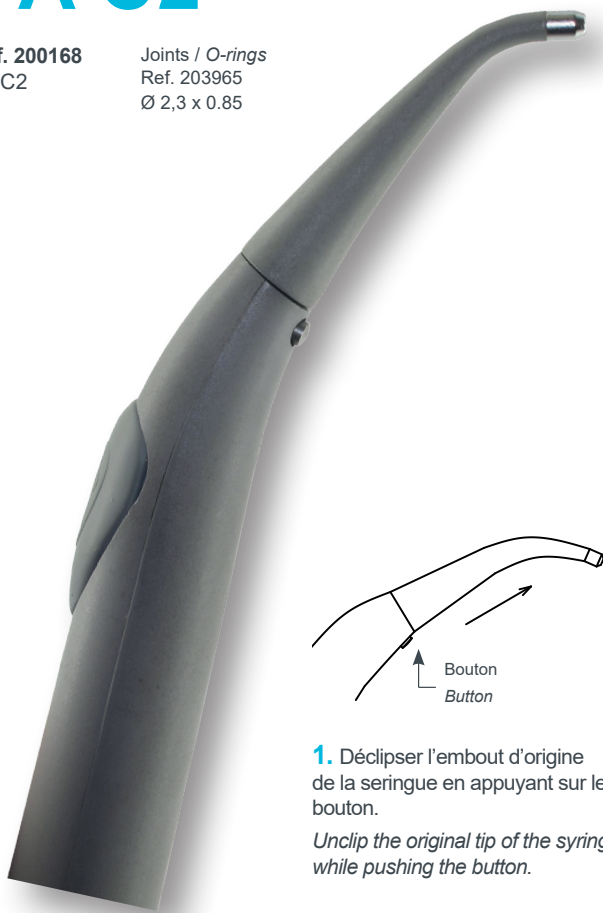
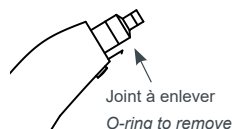
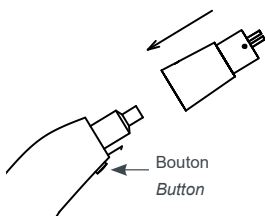
Unclip the original tip of the syringe while pushing the button.

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Remove the small O-ring situated at the end of the syringe and replace it by the one provided with the adapter.

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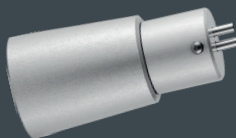
Lubricate the syringe, press the button of the syringe and push the adapter on the syringe while turning it.

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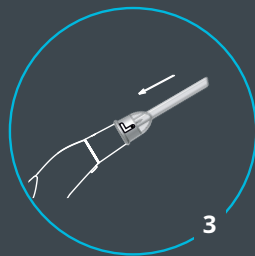
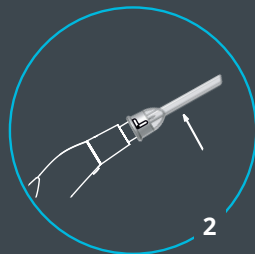
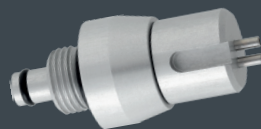
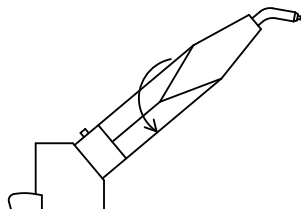
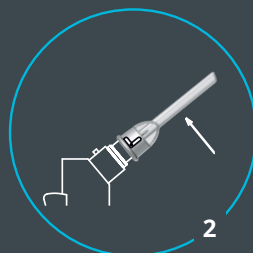
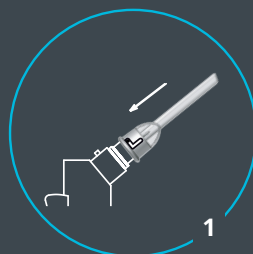
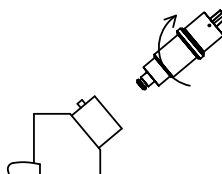
Lubricate the O-ring and screw the adapter.

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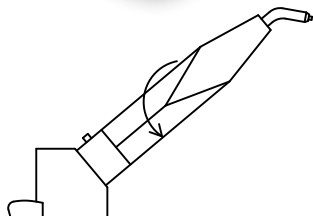
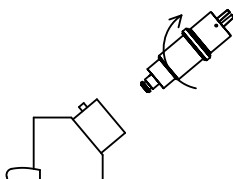
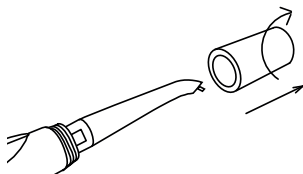
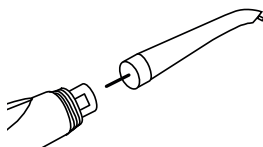
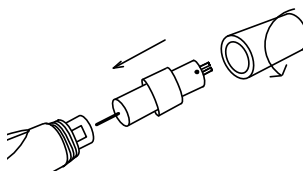
Unscrew the original nut from the syringe.

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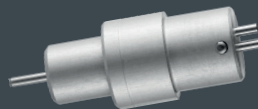
Pull the remaining part to remove it.

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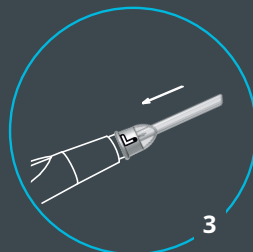
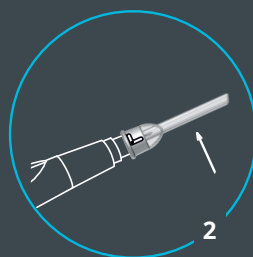
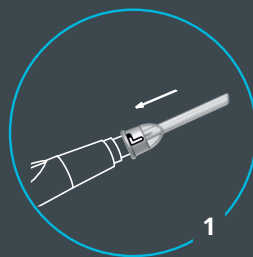
Push the adapter into the syringe, then screw the original nut of the syringe.

Adaptateur sans joint
Adapter without O-ring

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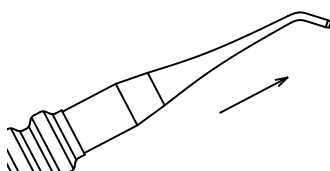
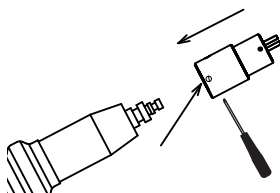
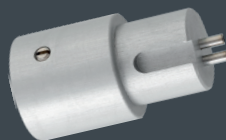
Unclip softly the original tip of the syringe. Lubricate abundantly the syringe.

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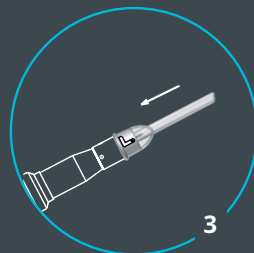
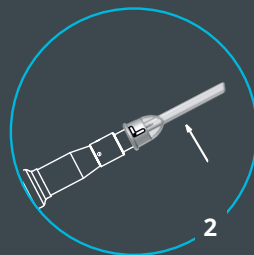
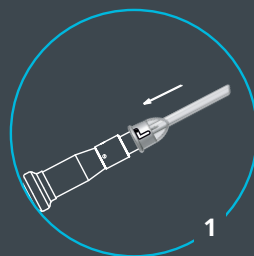
Unscrew the small screw then clip the adapter on the syringe while turning it and screw again the small screw.

Adaptateur sans joint
Adapter without O-ring

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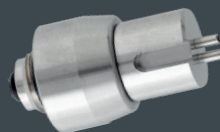
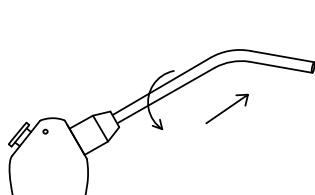
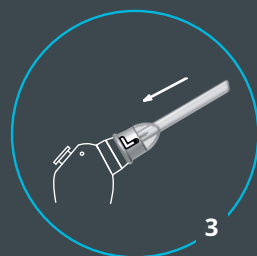
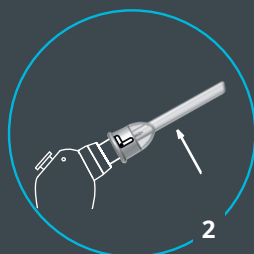
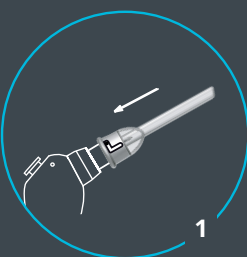
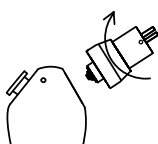
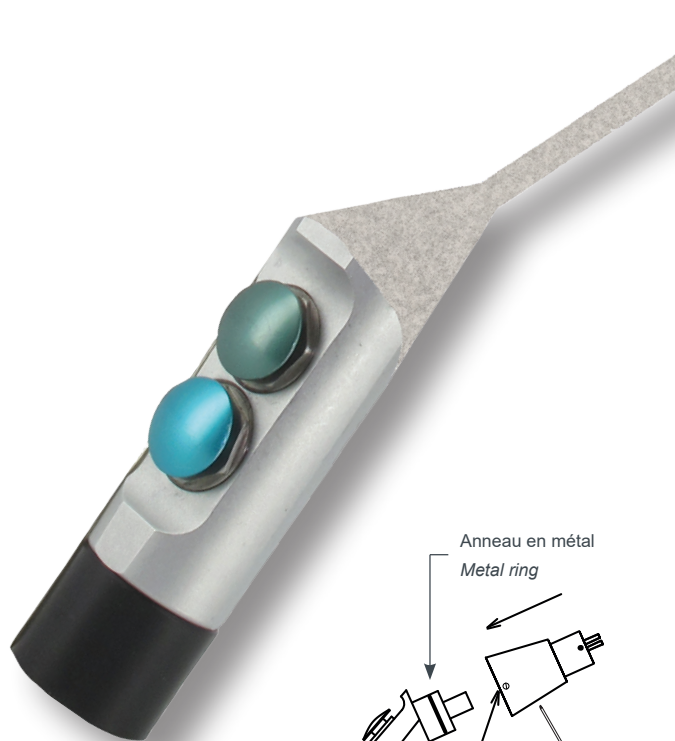
Lubricate the small O-ring at the end of the adapter then screw the adapter.

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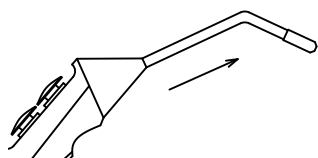
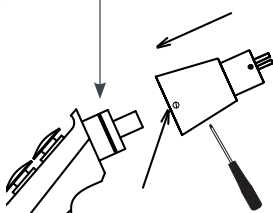
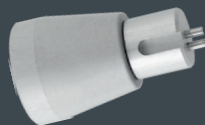
Unclip the original tip of the syringe.

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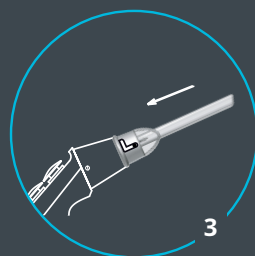
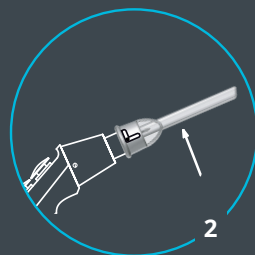
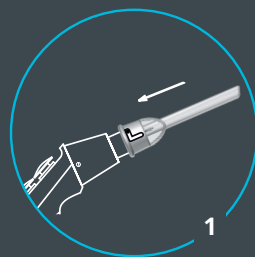
Check that the metal ring (the clip) is undamaged. Lubricate the syringe, unscrew the small screw, clip the adapter on the syringe while turning it, and screw again the small screw.

Adaptateur sans joint
Adapter without O-ring

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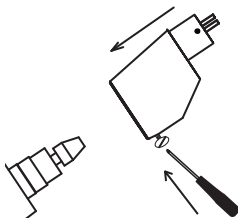
Unclip the original tip of the syringe. Check that the O-rings of the syringe are undamaged.

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Remove the small O-ring located at the end of the syringe (see drawing opposite) and replace it by the one provided with the adapter.

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Lubricate the syringe, clip the adapter on the syringe while turning it, then screw the little screw of the adapter.

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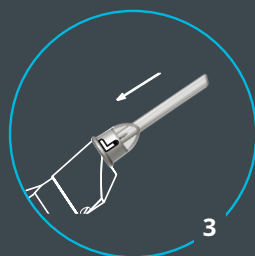
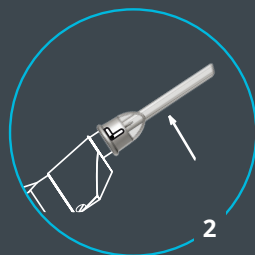
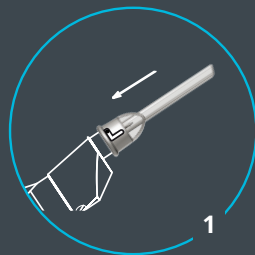
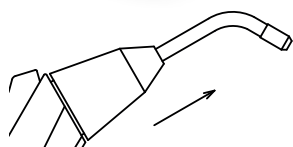
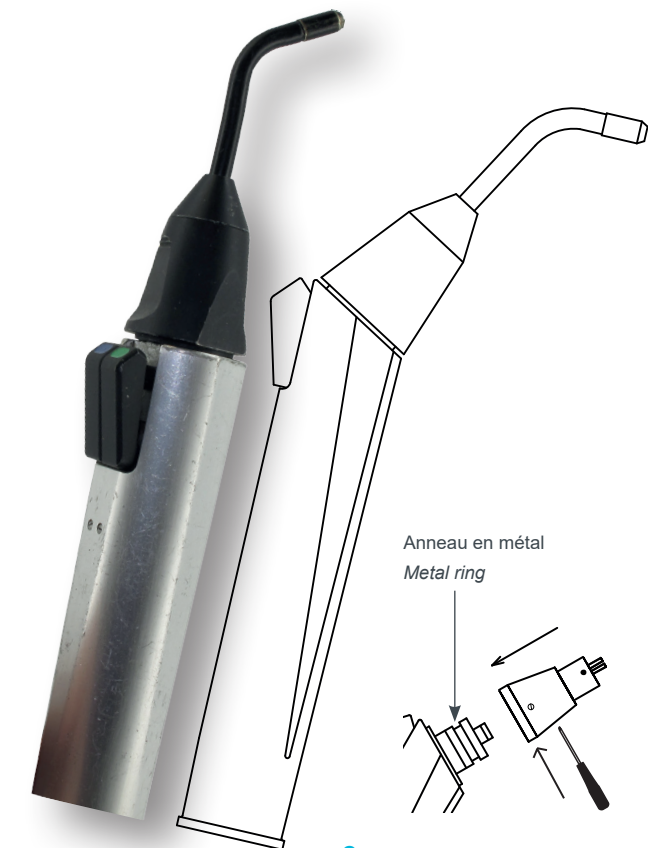
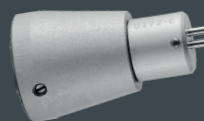
Check that the metal ring (the clip) is undamaged.

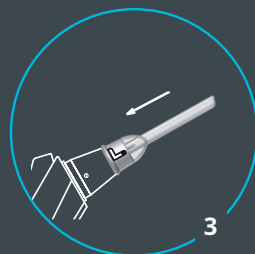
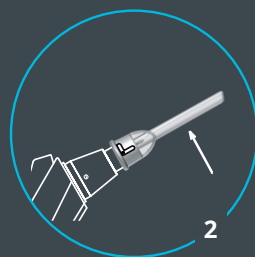
Lubricate the syringe, unscrew the small screw, clip the adapter on the syringe while turning it, and screw again the small screw.

Adaptateur sans joint
Adapter without O-ring

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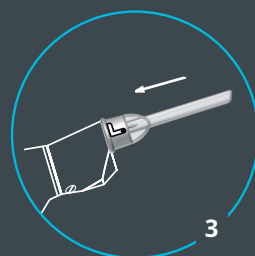
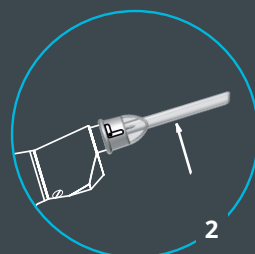
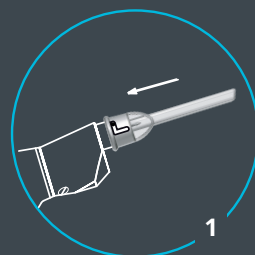
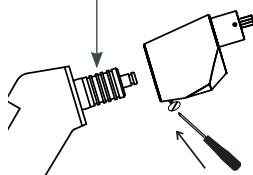
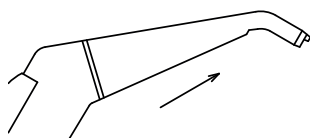
Check that the O-rings of the syringe are undamaged. Lubricate the syringe, unscrew the small crew, clip the adapter on the syringe while turning it and screw again the smal crew.

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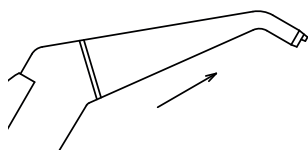
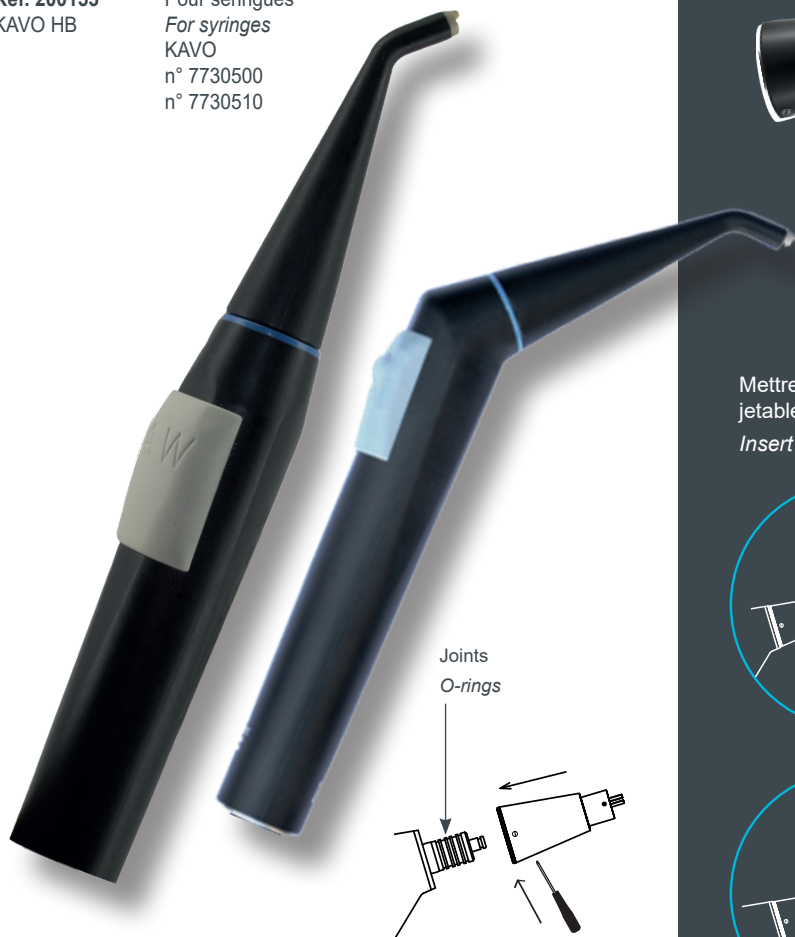
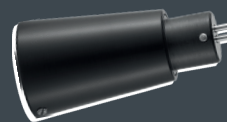
Check that the O-rings of the syringe are undamaged and lubricate them.

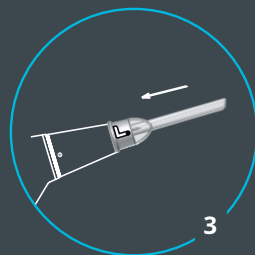
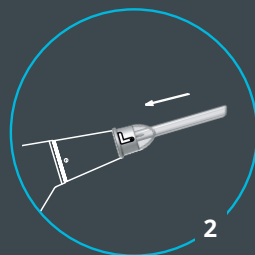
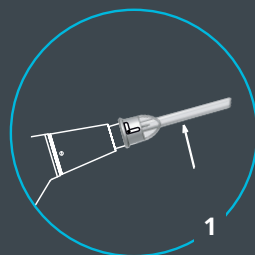
Unscrew the small crew, clip the adapter on the syringe while turning it and screw again the small crew.

Adaptateur sans joint
Adapter without O-ring

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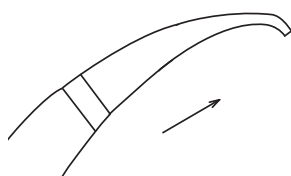
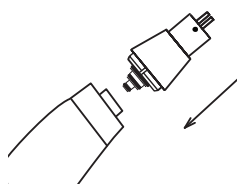
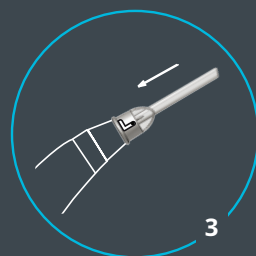
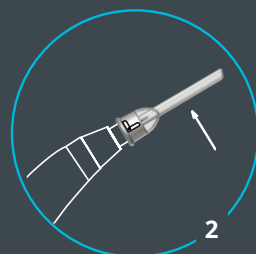
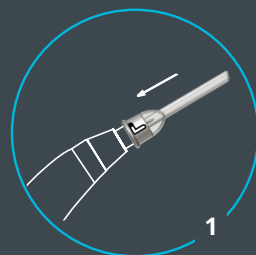
Lubricate the O-rings of the adapter. Clip the adapter on the syringe while turning it.

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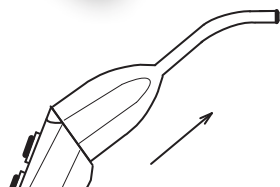
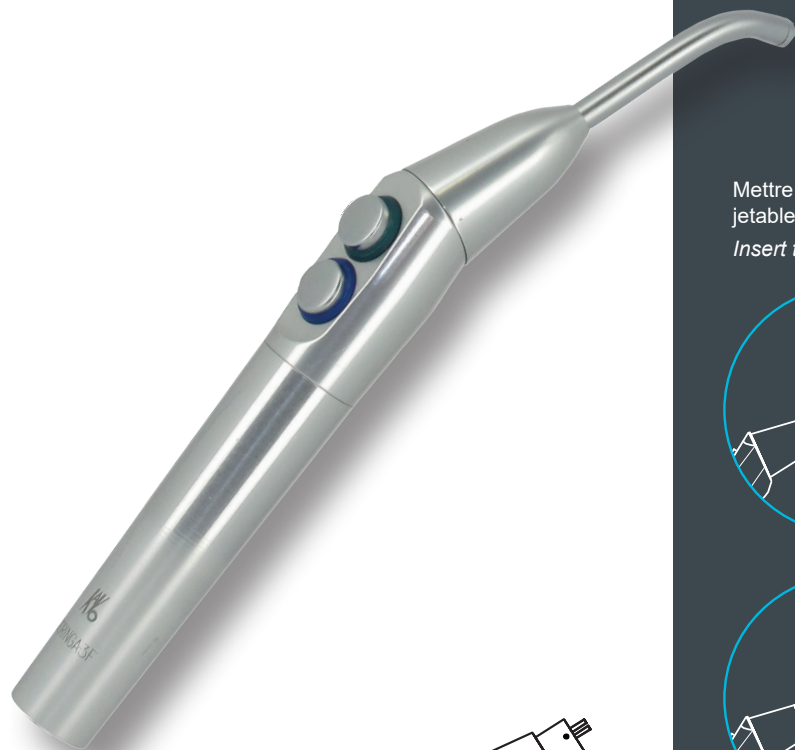
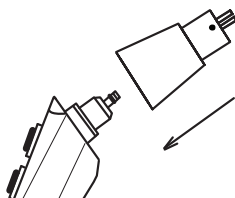
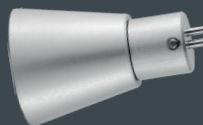
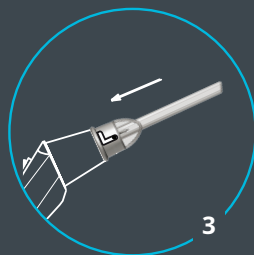
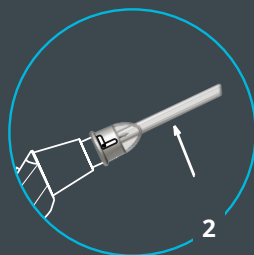
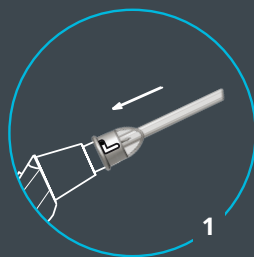
Lubricate the O-ring of the adapter. Clip the adapter on the syringe while turning it.

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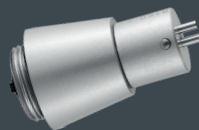
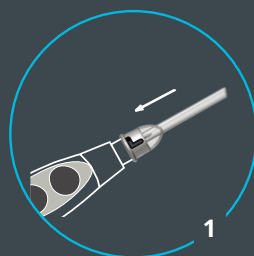
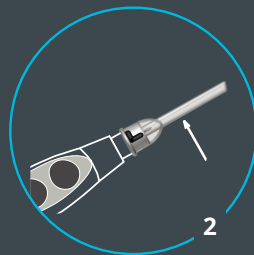
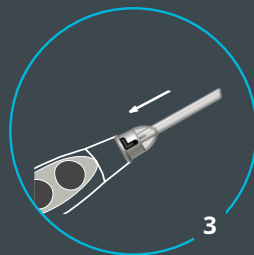
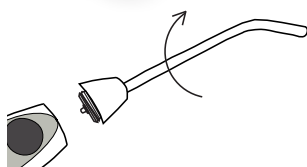
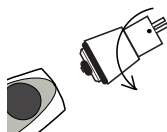
Lubricate the 2 O-rings of the adapter. Screw the adapter on the syringe while turning it.

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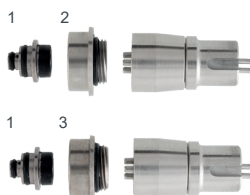
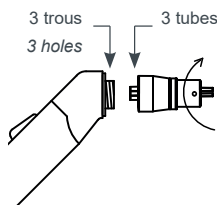
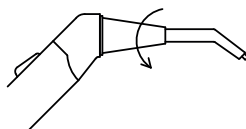
Warning: If your syringe is not anymore equipped with the original tip, remove the adapter already set on the syringe. Push (1) part 344 C inside the syringe, then screw (2) part 447 or part 222 (MINIMATE). Set the adapter LU in position as indicated above.

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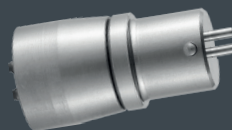
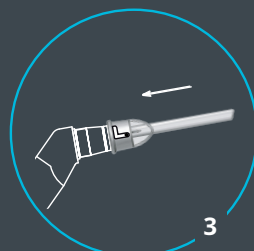
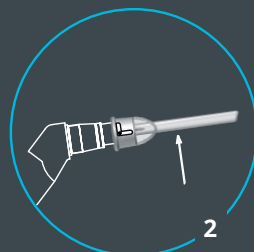
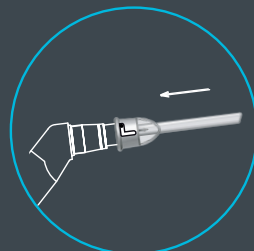
Put the 3 tubes of the adapter into the 3 holes of the syringe and screw the adapter.

Adaptateur sans joint
Adapter without O-ring

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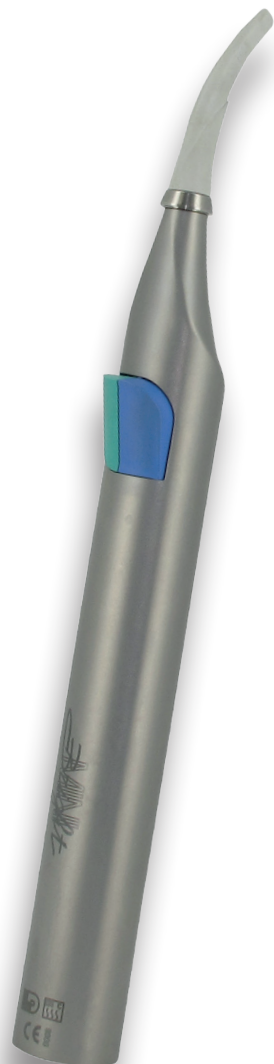
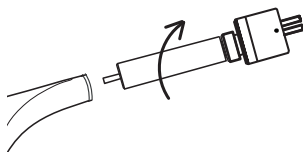
Unscrew the original tip of the syringe at the metal ring level.

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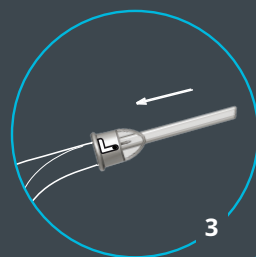
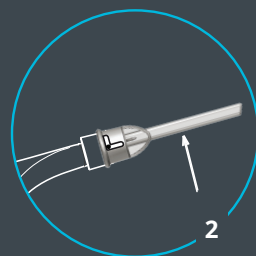
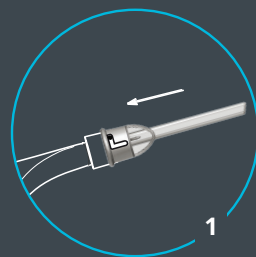
Insert all the way the adapter inside the syringe then screw it..

Adaptateur sans joint
Adapter without O-ring

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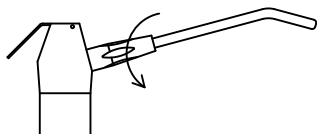
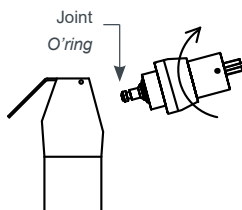
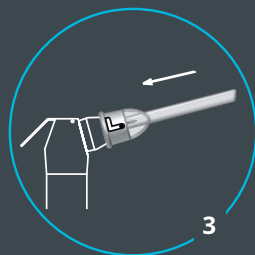
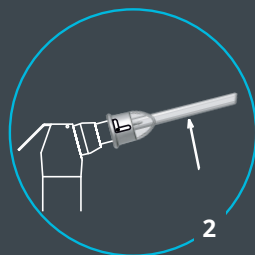
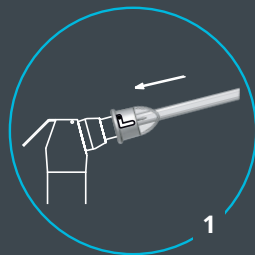
Check that the small O-ring is at the end of the adapter. Lubricate the O-rings and screw the adapter.

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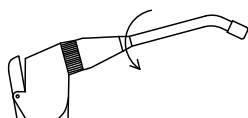
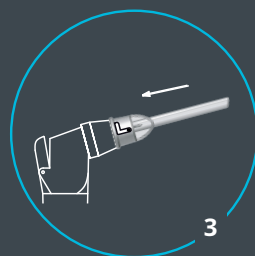
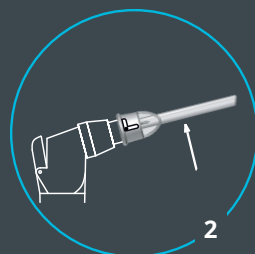
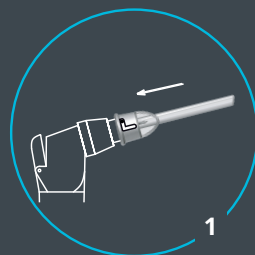
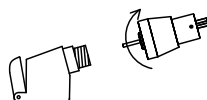
Lubricate the 2 O-rings and the part of the adapter going inside the syringe and screw the adapter.

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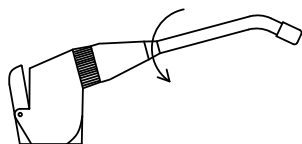
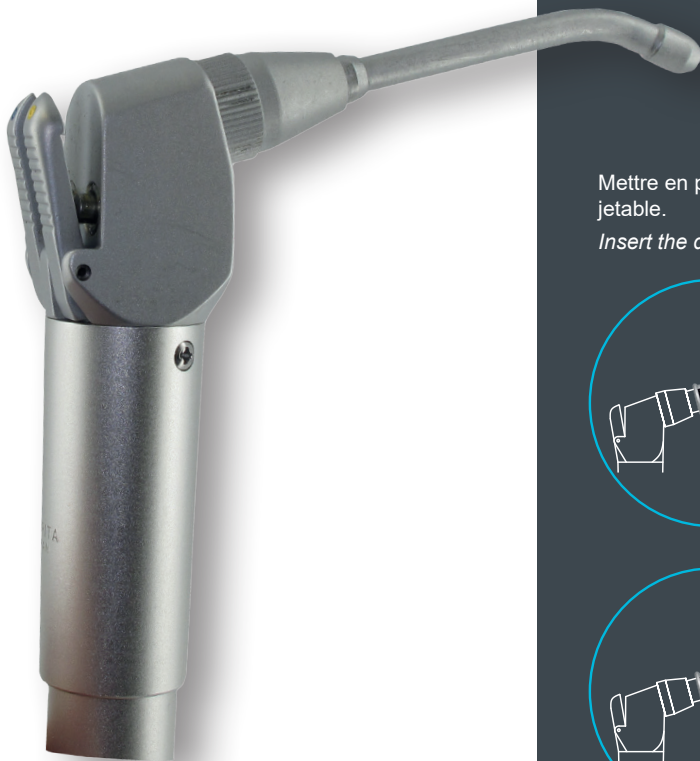
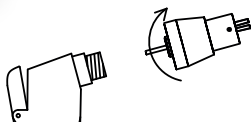
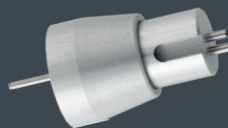
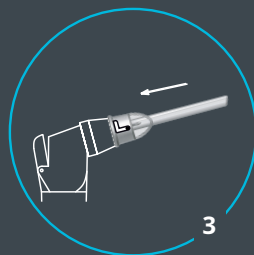
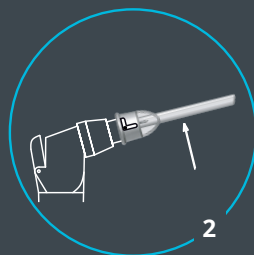
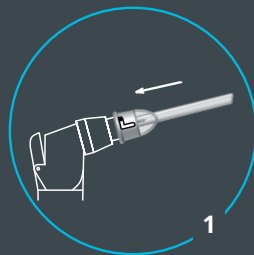
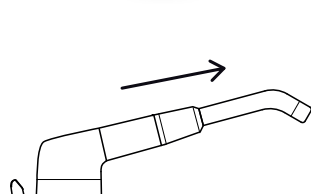
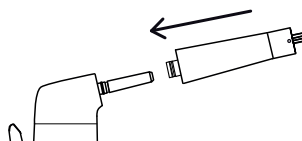
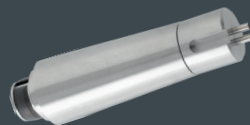
Unclip the original tip of the syringe

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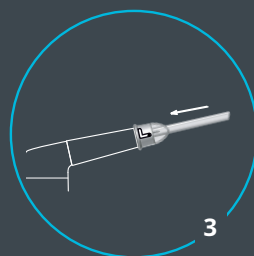
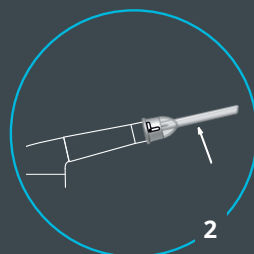
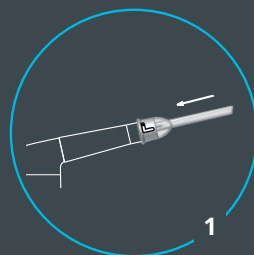
Lubricate the O-ring of the syringe and clip the adapter.

Adaptateur sans joint
Adapter without O-ring

Mettre en place l'embout jetable.

Insert the disposable tip.

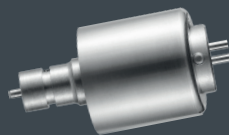
MORITA INOX WS97

Ref. 200174

MORITA INOX WS97

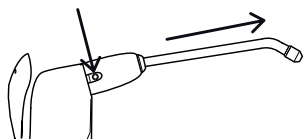
Adaptateur sans joint
Adapter without O-ring

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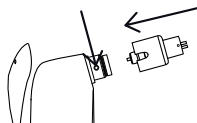
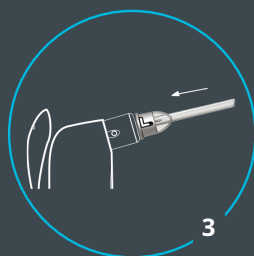
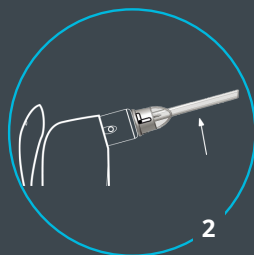
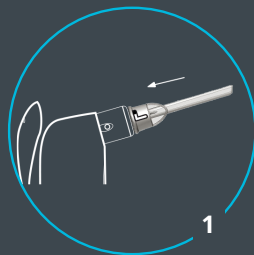
Press the button on the syringe and remove the original tip by pulling it.

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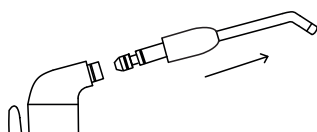
Press the button on the syringe and set up the adapter as shown in the drawing.

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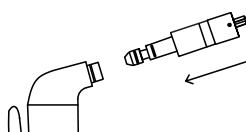
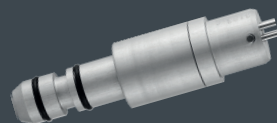
Pull on the original tip of the syringe to unclip it.

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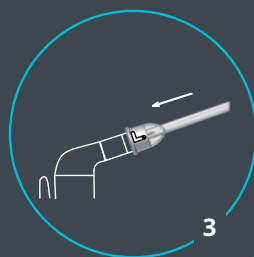
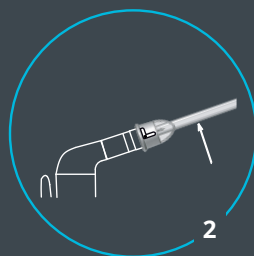
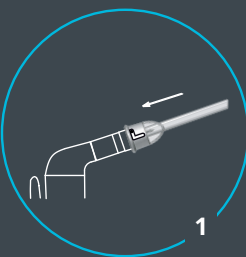
Lubricate the 2 O-rings of the adapter then clip it on the syringe.

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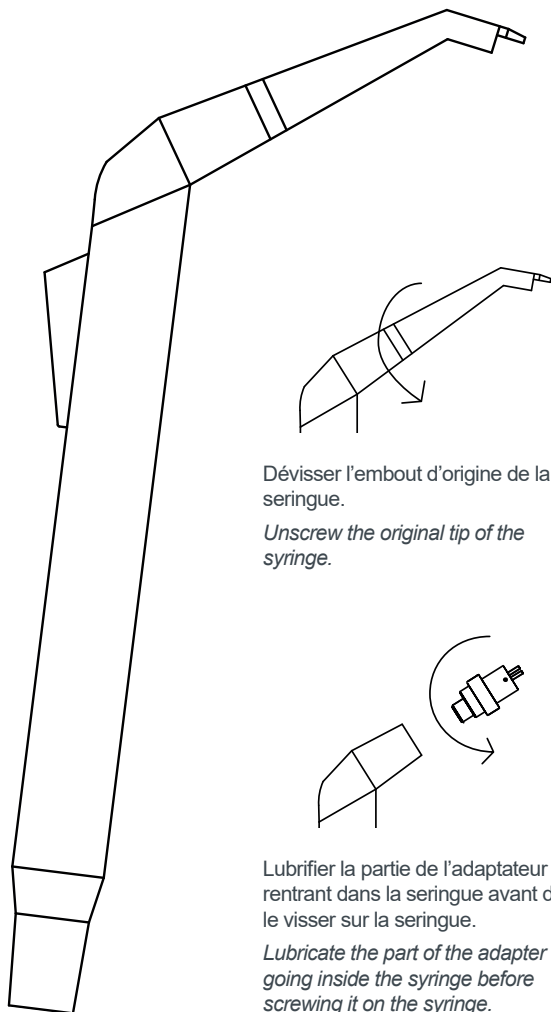
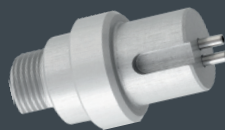
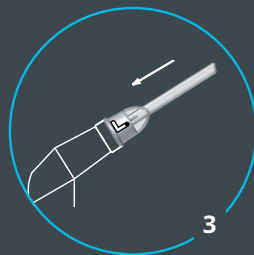
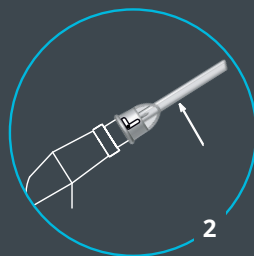
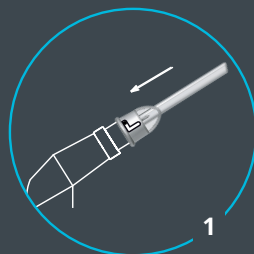
Lubricate the part of the adapter going inside the syringe before screwing it on the syringe.

Adaptateur sans joint
Adapter without O-ring

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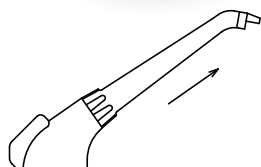
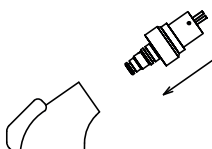
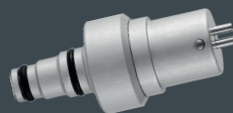
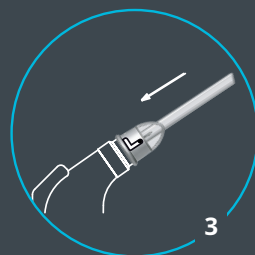
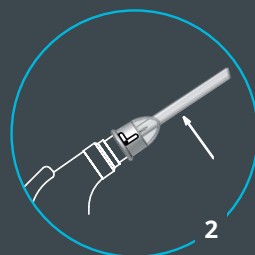
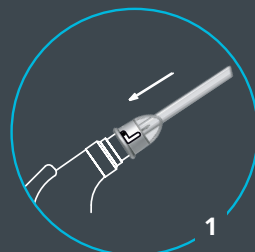
Lubricate the 2 O-rings and clip the adapter while turning it.

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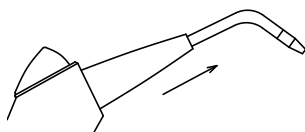
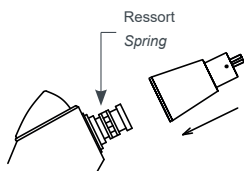
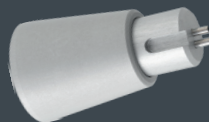
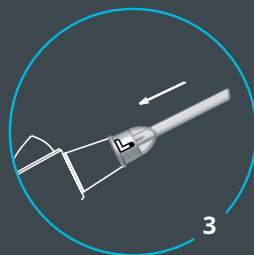
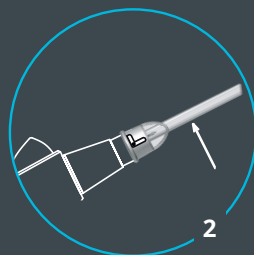
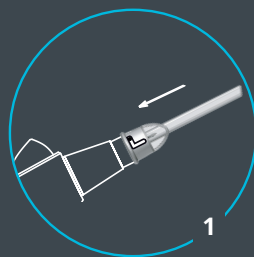
Check that the metal spring (the clip) is undamaged. Lubricate the syringe and clip the adapter while turning it.

Adaptateur sans joint
Adapter without O-ring

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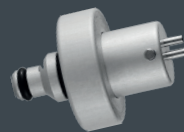
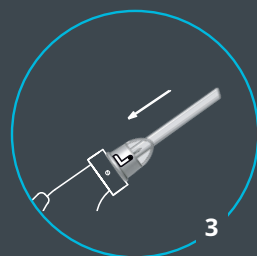
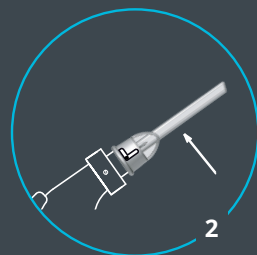
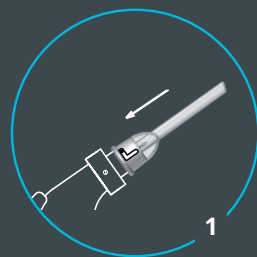
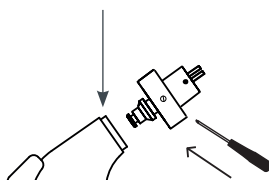
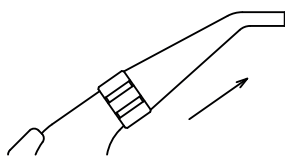
Check that the metal spring (the clip) is undamaged.

Unscrew the small screw, lubricate the 2 O-rings, push the adapter in the syringe while turning until you hear a clic and screw again the small screw.

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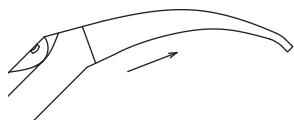
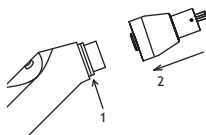
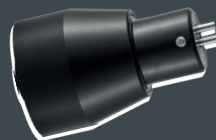
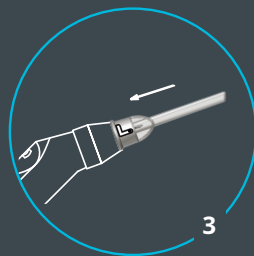
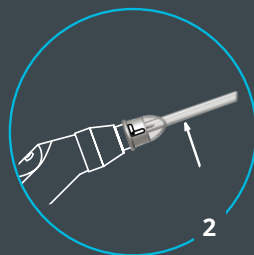
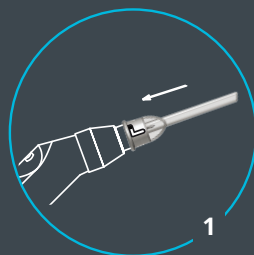
1. Check the good shape of the metal ring (1).
2. Lubricate the 3 O-rings and clip the adapter while turning (2). If you have a problem to clip it, you must change the metal ring of the syringe.

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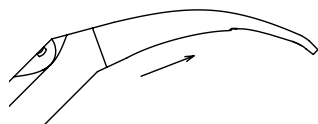
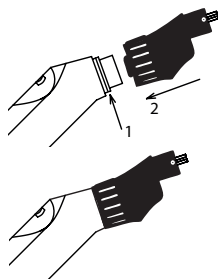
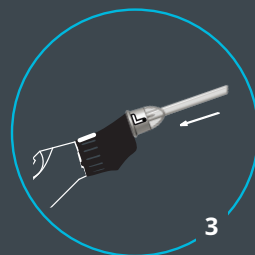
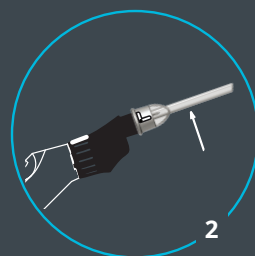
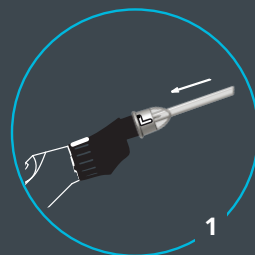
Check the good condition of the metal ring (1). Lubricate the 3 O-rings and clip the adapter while turning (2).

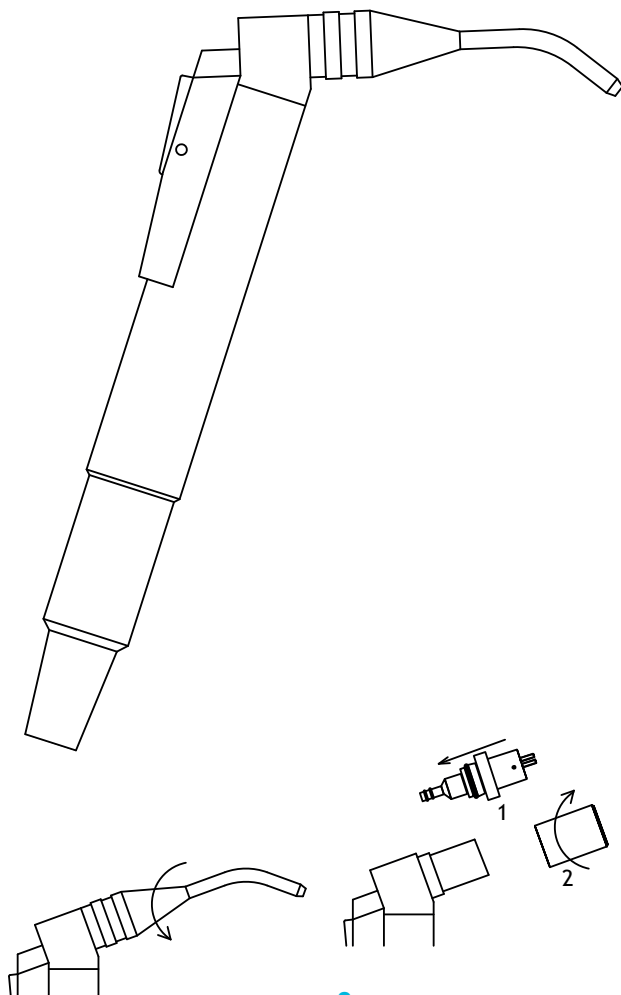
If you have a problem to clip it, you must change the metal ring of the syringe.

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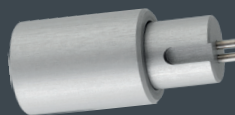
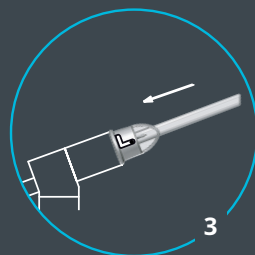
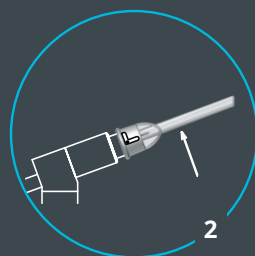
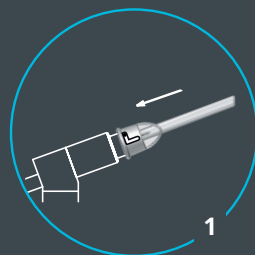
Lubricate the O-rings from the part n°1 and push it inside the syringe, then screw part n°2.

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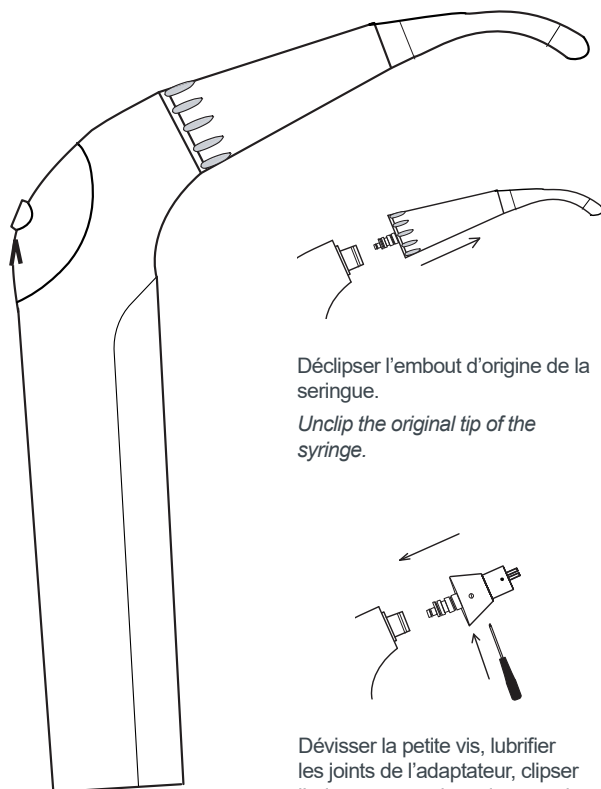
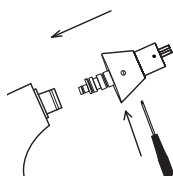
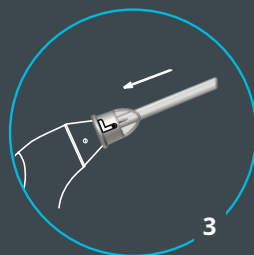
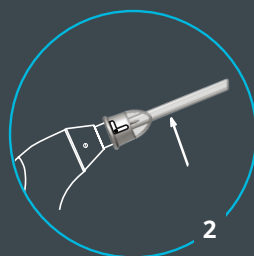
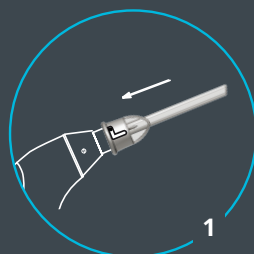
Unscrew the small screw, lubricate the O-rings of the adapter, clip the adapter on the syringe and then screw the small screw.

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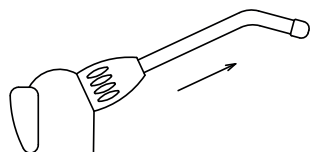
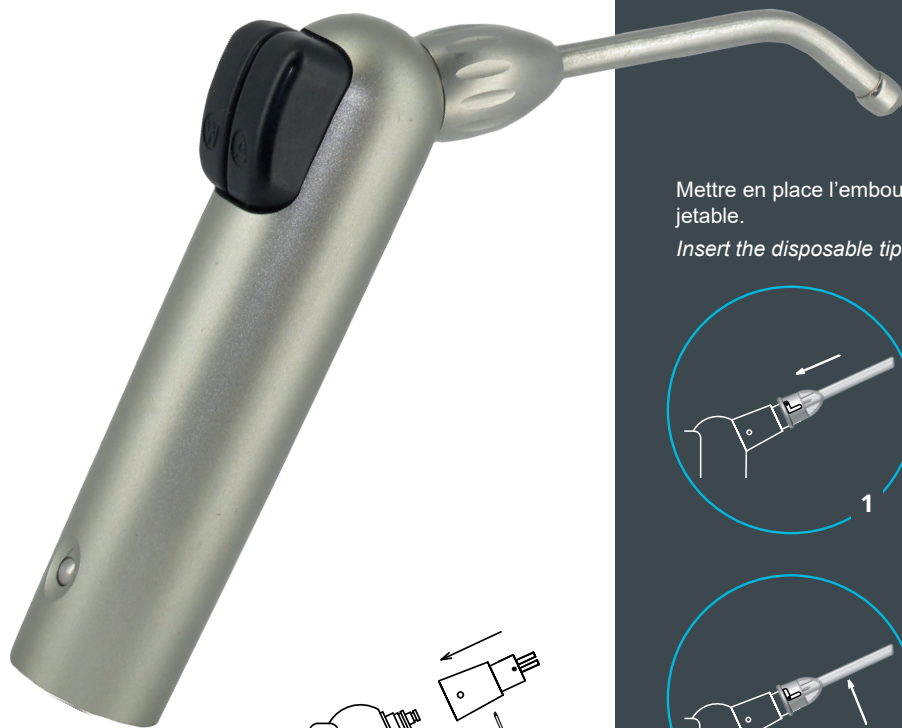
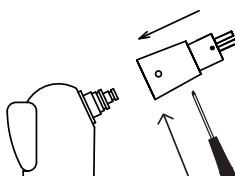
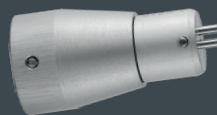
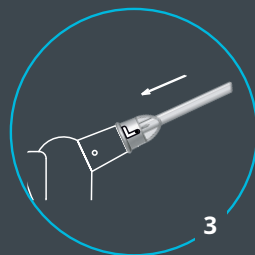
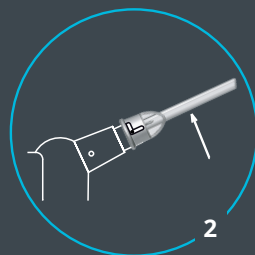
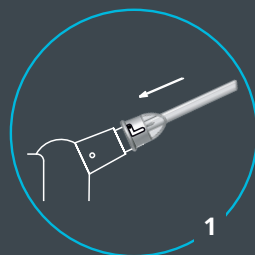
Lubricate the syringe, unscrew the small screw of the adapter then clip the adapter on the syringe while turning it and screw again the small screw.

Adaptateur sans joint
Adapter without O-ring

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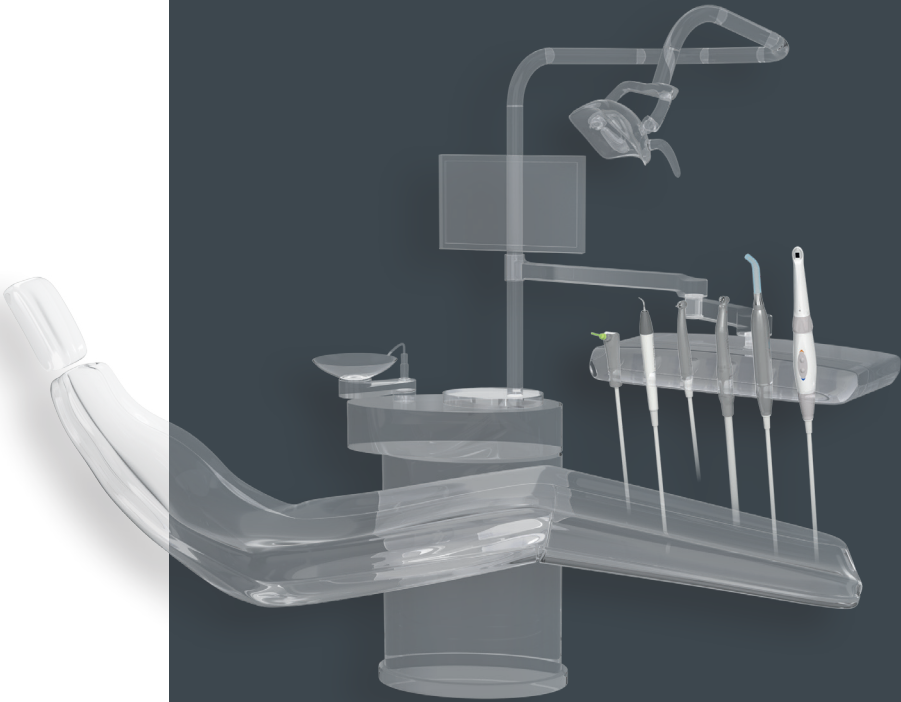
***Units, syringes brands and
corresponding adapters***

Fauteuil <i>Unit</i>	Modèle <i>Type</i>	Seringue <i>Syringe</i>	Adaptateur <i>Adapter</i>	Réf. <i>Ref.</i>
ADEC		Adec	Adec Inox	200195
A.D.M. Select Dental		Faro*	FA C ou/or FA C2	200167 200168
AIRCAR		Luzzani*	LU Inox	200150
AIREL		Faro*	FA C ou/or FA C2 ou/or FA S1	200167 200168 200182
		Luzzani*	LU Inox	200150
		Mini Bright	LU New Inox	200199
AJAX		Copie/Copy DCI	DCI	200152
		Copie/Copy Adec	Adec Inox	200195
AMPCO DENTAL		DCI	DCI	200152
ANATOM		Ritter	Ritter Topjet	200176
		Luzzani*	LU Inox	200150
ANCAR		Luzzani*	LU Inox	200150
ANTHOS		Luzzani*	LU Inox	200150
		Faro*	FA S1	200182
		Faro*	FA E	200165
		Anthos	Cefla Inox	200197
		Nlle/New Cefla	Cas Inox	200198
ASISTROM		Luzzani*	LU Inox	200150
AS GOLD		DCI	DCI	200152
ASEPTICO		Adec	Adec Inox	200195
		Forest	Adec Inox	200195
B-PRO- DUCTION		Copie/Copy DCI	DCI	200152
BEAVER DENTAL		Adec	Adec Inox	200195

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*Faro and Luzzani are not trademarks of ACTEON Manufacturing (ACTEON® Group)

Fauteuil Unit	Modèle Type	Seringue Syringe	Adaptateur Adapter	Réf. Ref.
BELMONT		Belmont	Belmont New	200158
		Faro*	FA C ou/or FA C2	200167 200168
		Adec	Adec Inox	200195
		Luzzani*	LU Inox	200150
		DCI	DCI	200152
BIO-DENTAL		Faro*	FA C ou/or FA C2	200167 200168
BRAVO 21C		DCI	DCI	200152
BROX-DENT		DCI	DCI	200152
CACCAN		Forest	Adec Inox	200195
		Faro*	FA C ou/or FA C2	200167 200168
		Luzzani*	LU Inox	200150
CANCAN		DCI	DCI	200152
CASMON		Copie/Copy DCI	DCI	200152
CASTEL-LINI	SKEMA CLEAN WARM 2	Castellini	Castellini V	200163
	LOGOS	Castellini	Castellini Logo New	200184
	PUMA, PUMA+ Threesteril LOGOS Junior	Castellini Logo (boutons ronds / round buttons)	Castellini Logo New	200184
	LOGOS Junior, Duo	Avec lumière <i>With light</i>	Cefla Inox	200197
		Nlle/New Cefla	Cas Inox	200198

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INCIDENT	DETECTION	SOLUTION
<i>Difficulty to set up the adapter</i>	<i>Check you have the right adapter to fit your type of syringe. Check the correct set up of the adapter on the syringe.</i>	<i>Contact your ACTEON® salesperson. Read carefully the assembly instructions.</i>
<i>Difficulty to fit up the adapter on the syringe</i>	<i>Check the O-rings of the adapter and the syringe have been lubricated. Check the good shape of the clips (metal rings) of the syringe.</i>	<i>Lubricate the O-rings before setting up the adapter. Ask the person in charge of the maintenance of your unit to change this part.</i>
<i>Loss of power of the air-water jet after setting up the adapter</i>	<i>Check the air and water pressure of the unit (the difference between the water and air pressure should not exceed 1 bar).</i>	<i>Set the air pressure of the unit at 3.5 bars minimum (about 50 PSI).</i>
<i>Water through the air jet of the syringe</i>	<i>Check the good setting of the water and air pressure of the unit (the difference between the water and air pressure should not exceed 1 bar).</i>	<i>The air pressure of the unit must be set at 3.5 bars minimum (about 50 PSI).</i>
<i>Water leak between the adapter and the disposable tip</i>	<i>Check the disposable tip is set correctly on the adapter.</i>	<i>Set correctly the tip on the adapter.</i>
<i>Bad spray, not enough strong</i>	<i>Check the good setting of the water and air pressure of the unit (the difference between the water and air pressure should not exceed 1 bar).</i>	<i>The air pressure of the unit must be set at 3.5 bars minimum (about 50 PSI).</i>
<i>Absence of water or air (or both) discharge from the tip</i>	<i>Check the setting up of the tip on the adapter.</i>	<i>Clear the 2 tubes of the adapter with an instrument (barbed broach).</i>

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INCIDENT	DETECTION	SOLUTION
<i>Problem of water in the air after the set up of the adapter</i>	<i>Check the good shape of the O-rings of your syringe.</i>	<i>Change the O-rings of the syringe in bad shape.</i>
<i>Problem of water in the air tube after setting up of the adapter ADEC Inox or any other adapter screwed on the syringe</i>	<i>Check the positioning and the good shape of the small O-ring at the end of the adapter.</i>	<i>Unscrew the adapter, set the small O-ring at the right place, change it if it is in bad shape and screw again the adapter.</i>
<i>Problem of water in the air tube after setting up of the adapter for DCI syringe</i>	<i>Check the good shape of the O-ring inside the syringe.</i>	<i>Change the O-ring inside the syringe.</i>
<i>Problem of water in the air on the syringe KAVO DB, Kavo GB, Kavo HB</i>	<i>Check the good shape of the O-rings.</i>	<i>Change the O-rings of the syringe in bad shape.</i>
<i>Sirona B onto the latest Sirona syringe (black with blue buttons)</i>	<i>Verify the good shape of the metal ring on the syringe (the clip).</i>	<i>Change the metal ring of the syringe if the shape is not good.</i>

**ACTEON
Manufacturing**

A company of ACTEON® Group

17 avenue Gustave Eiffel,
ZI du Phare,
33700 Mérignac, France

Tel. +33 (0) 556 340 607
Fax. +33 (0) 556 349 292

info@acteongroup.com

Follow-us