


Main unit 1 unit	65um supragingival sodium bicarbonate powder 1 bottle
Foot pedal 1 pc	25um subgingival glycine powder 1 bottle
Power adapter (power cord included) 1 pc	Subgingival nozzle 10 pcs
Ultrasonic handpiece (HW-5L) 1 pc	Nozzle wrench 1 pc
Supragingival air polishing handpiece (AP-1) 1 pc	Sterilization box 2 pcs
Subgingival air polishing handpiece (AP-2) 1 pc	1400mL water bottle 1 pc
Tips P1×2 P3×1 P8×1 G1×1 G4×1	700mL conical water bottle 1 pc
1# Torque wrench 1 pc	Needles (Long & Short) 1 kit
Supragingival powder tank 1 pc	Union tee 1 pc
Subgingival powder tank 1 pc	Male and female connectors 1 pair
Air pipe connector 1 pc	

Guilin Woodpecker Medical Instrument Co.,Ltd.
Information Industrial Park,Guilin National High-Tech
Zone,Guilin,Guangxi, 541004 P.R.China
Sales Dept.: +86-773-5873196
E-mail: woodpecker@glwoodpecker.com
Website: http://www.glwoodpecker.com
ZHEWVC-072

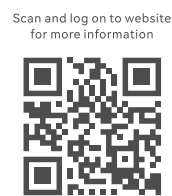

AP-B Dental Scaler and Air Polisher

Faster scaling

Faster and more convenient scaling



Optimized piping structure

New aerodynamic design
Reduce fluid pressure loss
Effectively remove plaque and pigment



Before scaling



Plaque disclosing



After scaling

Classic sodium bicarbonate powder

Average particle size:

about 65 microns; hardness is less than that of enamel, doing no harm to teeth

Contain natural plant-derived essence for a more pleasant taste

The powder is waterproof to prevent the device from being blocked by wet powder

Teeth	Moh's hardness
Enamel	4-5
Dentin	2-2.5
Cementum	2-4

Powder	Moh's hardness
Sodium bicarbonate	2.5
Glycine	2
Calcium carbonate	3
Aluminum hydroxide	3

The promotion and sale of SUPER and SUPER PRO are not aimed at the following countries: France, Spain, Germany, Switzerland, Norway, Czech Republic, United Kingdom, Italy, Netherlands, Austria, Belgium, Denmark, Turkey, Sweden, Liechtenstein, Japan, and the United States.



AP-B

Stand out
from the crowd

Innovative design

Soft candy color

Stand out
where others
fade into the masses

A unique Dental Scaler
and Air Polisher is in!



High-performance ultrasonic scaling

High-quality and efficient tooth cleaning adjunct
Easy removal of large pieces of calculus
Detachable and autoclavable ultrasonic handpiece


Automatic identification of scaling mode

Pick up the Air Polishing Handpiece to automatically enter the Air Polishing Mode

Pick up the Ultrasonic Handpiece to automatically enter the Ultrasonic Scaling Mode




Air
Polishing
Mode

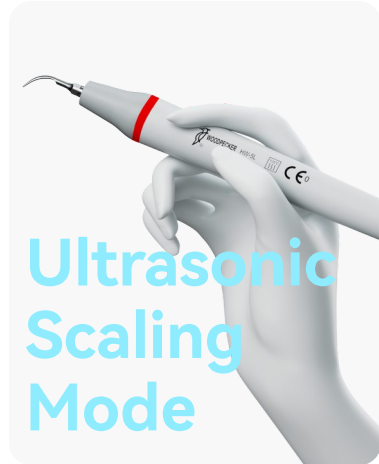

Ultrasonic
Scaling
Mode


Two-in-one: Supragingival & subgingival mode

Three-section handpiece design

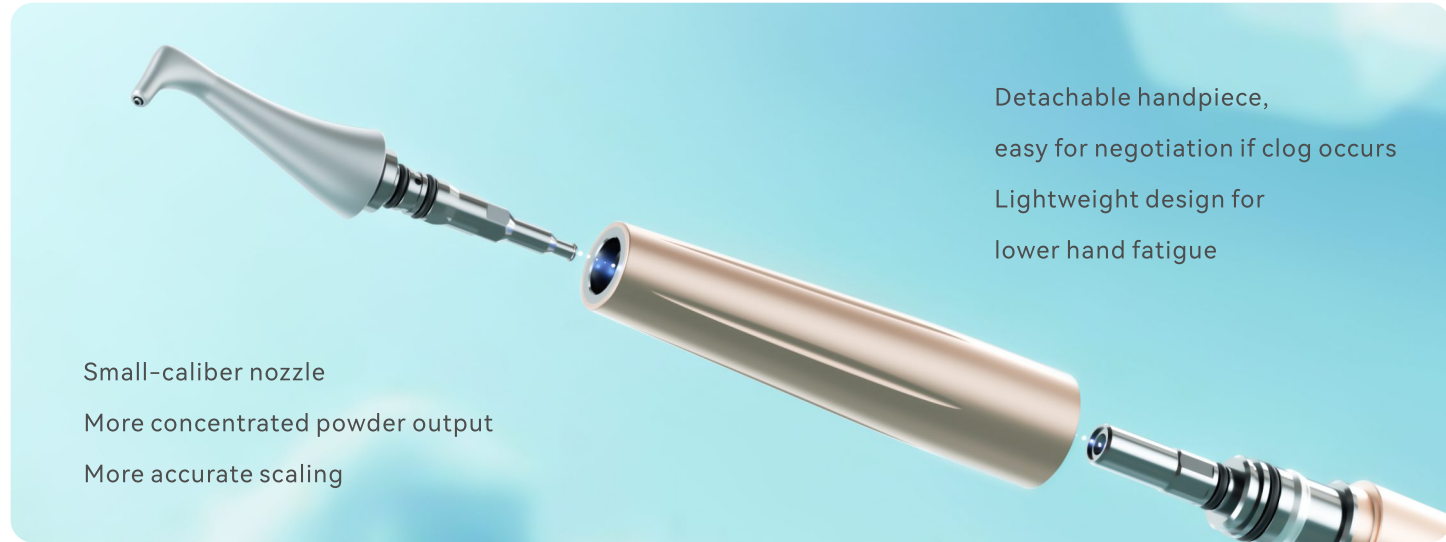

Small-caliber nozzle
More concentrated powder output
More accurate scaling

Detachable handpiece,
easy for negotiation if clog occurs
Lightweight design for
lower hand fatigue

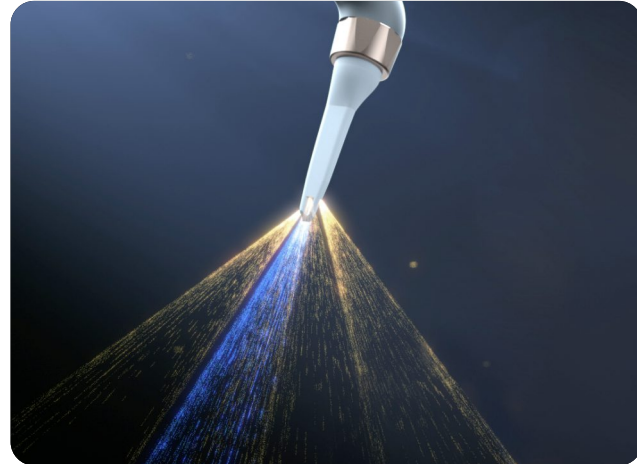

Subgingival nozzle made of medical polymers
360° rotatable
Bendable

Three-direction powder output
One-direction water output
Thorough plaque removal
Great accessibility

Ergonomic design

Fully transparent funnel-shaped powder tank

Stable powder output

360°rotatable, allowing real-time observation of powder flow in the tank

Compact and lightweight

The main unit weighs only 2kg

Space-saving base design

