




SOOO ADHERENT  
IS **MUCOPREN**<sup>®</sup> SOFT

SOFT RELINING

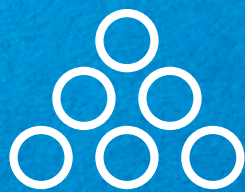

**KETTENBACHDENTAL**

Simply intelligent



# MUCOPREN® SOFT – LIVES UP TO ITS PROMISE

## Denture relining

Mucopren® Soft is a permanently soft durable reline long-term material based on A-silicone for direct and indirect application.



## INDICATIONS

- for permanently soft linings
- for compensating for areas of mucous membrane with varying degrees of resiliency
- for relieving sharp bony ridges
- for relieving the operation site after placing implants
- for relining temporary restorations
- for conditioning the mucous membranes in cases with denture sores
- for fabricating and relining obturators

## Technical data

Product name	Mixing ratio	Working time at 23 °C / 74 °F ≤	Working time at 35 °C / 95 °F ≤	Intraoral setting time at 35 °C / 95 °F ≥	Total setting time ≥	Hardness (approx.) Shore	Linear dimensional change (approx.) %	Elastic recovery test (approx.) %	Strain in compression (approx.) %
Mucopren® Soft direct relining	1:1	2:15 min.	–	3:15 min.	5:30 min.*	A 28	– 0.20	99.5	9.4
Mucopren® Soft indirect relining	1:1	2:15 min.	–	–	> 30 min. in pressure vessel	A 28	– 0.20	99.5	9.4

\*Total setting time (removal time from mouth) from start of mixing

## INHIBITS BACTERIAL COLONIZATION



- Thanks to its smooth, hydrophobic silicone surface\*

## EASY PROCESSING



- Due to good millability

## HIGH WEARING COMFORT

- As a result of permanently elastic properties

## GUARANTEED LONG-TERM ADHESION TO THE DENTURE

- Thanks to the excellent adhesion of the **Mucopren®** Adhesive





# MUCOPREN® SOFT – THE RANGE

## Mucopren® Soft

### Base set

50 mL relining silicone, 50 mL silicone sealant,  
10 mL adhesive, 7 green mixing tips, 20 blue mixing tips, 1 steel bur,  
1 brush holder, 20 single-use brushes, accessories

### Relining silicone

2 × 50 mL cartridges, 6 mixing tips

### Mucopren® Silicone sealant

50 mL cartridge, 10 mixing tips

### Mucopren® Adhesive

10 mL bottle



Kettenbach GmbH & Co. KG  
Im Heerfeld 7  
35713 Eschenburg, Germany  
Phone +49 (0)2774 7050  
info@kettenbach.com  
www.kettenbach-dental.com



022746\_0222

  
**KETTENBACHDENTAL**  
Simply intelligent