



# kit and refills

## myQuickmat hoop kit

circumferential matrix system



REF 6600

- 40 pcs. Quickmat hoop, preloaded circumferential matrices, assorted (20 x regular, 20 x tabbed)
- 2 pcs. myHooper instrument
- 100 pcs. myHoopy wedge, interdental plastic wedges, assorted



for more information

### myHooper instrument

Contents	REF
1 pc. myHooper instrument	6602



### Quickmat hoop

preloaded circumferential matrices

Contents	Size	REF
50 pcs. Quickmat hoop, regular	small 5.5 mm	6603
50 pcs. Quickmat hoop, regular	small 6.5 mm	6604
50 pcs. Quickmat hoop, regular	large 5.5 mm	6605
50 pcs. Quickmat hoop, regular	large 6.5 mm	6606

Contents	Size	REF
45 pcs. Quickmat hoop, tabbed	small 5.5 mm	6607
45 pcs. Quickmat hoop, tabbed	small 6.5 mm	6608
45 pcs. Quickmat hoop, tabbed	large 5.5 mm	6609
45 pcs. Quickmat hoop, tabbed	large 6.5 mm	6610
45 pcs. Quickmat hoop, tabbed	23 small, 22 large 8.0 mm	6611

### myHoopy wedge

interdental plastic wedges

Contents	Size	REF
100 pcs. myHoopy wedge	yellow (S)	7602
100 pcs. myHoopy wedge	green (M)	7603
100 pcs. myHoopy wedge	teal (L)	7604
100 pcs. myHoopy wedge	blue (XL)	7605
100 pcs. myHoopy wedge, assorted	(30xS, 30xM, 20xL, 20xXL)	7601



MKTG-069-FL-EN REV. 06.2026

# myQuickmat hoop

## new! circumferential control you can count on



### ADVANTAGES

easy placement similar to a Tofflemire retainer – without its drawbacks

detachable instrument for fast placement, tightening and adjustment

anatomically contoured band for leak-free margins and natural restorations

compact design that is ideal for tight spaces, including second and third molars

35-micron hard contoured foil: strong, shape-retentive and passes tight contacts easily

colour-coded adjustable plastic retainer that stays secure without loosening

swivelling retainer arm for use in any quadrant

Polydentia SA  
via Cantonale 47  
6805 Mezzovico-Vira  
Switzerland

phone: +41 91 946 29 48  
info@polydentia.ch  
polydentia.ch



## circumferential matrix system





new!

# myQuickmat hoop

circumferential matrix system

circumferential control you can count on




"The Quickmat hoop matrix is extremely user-friendly: even extensive composite restorations can be completed quickly, with excellent predictability of the final result. Interproximal contacts and tooth anatomy are preserved, and the time required for the final finish is always minimised."

Dr. Thomas Casanova - Switzerland

The myQuickmat hoop kit is a circumferential matrix system that preserves the familiar placement workflow of a Tofflemire retainer, while eliminating many of its traditional constraints. With easy placement, anatomical contours, a compact design, hard foil, and full swivelling capability, myQuickmat hoop sets a new benchmark for modern circumferential matrix band systems.

The system includes a selection of regular and tabbed preloaded circumferential matrices, an instrument and interdental plastic wedges, which are useful for stabilising the matrices and separate the adjacent tooth.

Clinical image by courtesy of Dr. Perucchi & Dr. Scognamiglio, Switzerland

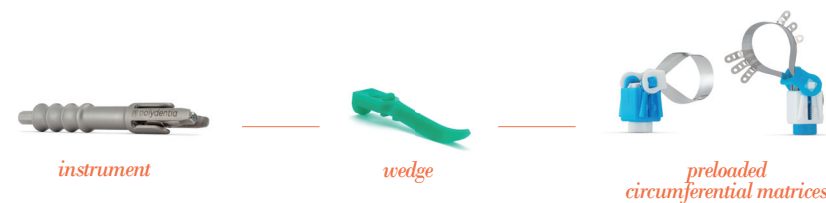

## Quickmat hoop

The Quickmat hoop preloaded matrix is a contoured circumferential matrix designed to deliver precision and consistent clinical results. Its detachable instrument (myHooper instrument) allows fast, intuitive placement using a Tofflemire-style technique but without leaving a bulky holder in the patient's mouth. This eliminates interference with soft tissues and prevents the matrix from being displaced during the procedure. The colour-coded adjustable retainer remains securely in place without loosening, while the swivelling arm allows use in any quadrant.

Quickmat hoop preloaded matrices are available regular and tabbed, in two circumferential sizes (small and large) and three height options (5.5, 6.5, 8.0mm).

The Quickmat hoop matrices are available in a tabbed version to help with large multi-surface composite restorations. The featured Control Hoop Tabs are specifically designed to facilitate reliable tight contacts when the matrix band unintentionally pulls away from adjacent teeth and open contacts occur.



regular





tabbed

The shape at the base of the band nips in and then flares out. It imitates the shape of a regular tooth which, obviously, is not straight up and down.



To maintain firm contact and ensure tight interproximal contacts:

- **option 1**  
Pull on the tab while curing the proximal composite.
- **option 2**  
Use flowable composite to tack the tabs to the adjacent tooth.

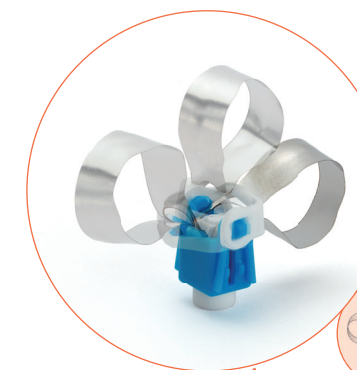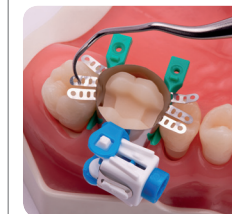

quadrant swivel

swivelling arm allows it to be used in any quadrant



myHooper instrument detachable instrument

myHoopy wedge

hard contoured 35µm Foil Band


colour-coded, adjustable retainer



compact design

## myHoopy wedge

The myHoopy wedge provides reliable separation, effective tissue compression and a secure matrix seal to help prevent seepage. Its parallel walls and gradually increasing rigidity allow for smooth insertion and firm stabilization once in position.

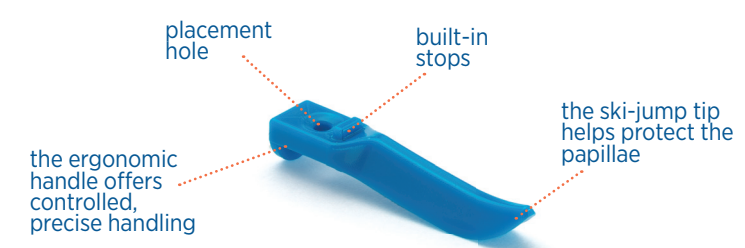

placement hole

built-in stops

the ergonomic handle offers controlled, precise handling

the ski-jump tip helps protect the papillae