

#whdentalwerk



video.wh.com

Now at
your dealer
or wh.com



Complete and versatile: Accessories



teon
teon +

Complete your device with the interior you really need

Baskets



A102 Upper basket, 19721400
Loading: H 205 (+/-30), W 475, D 443 mm



A105/1 Upper basket/injectors, 19721455
Loading: H 235, W 360, D 435 mm



A151 Lower basket, 19721402
Vertical clearance: H 230 (+/-30)/H 210 with A102/A105/1. Loading: H 88, W 529, D 522 mm

Connectivity



XKM RS232 10 Med Serial module, 19721432
Serial module



XKM 3000L Med Ethernet module, 19721433
Ethernet module

Adaptors



A801 Adaptor, 19721405
For use of AUF1 in combination with A105/1.



A803 Holder, 19721411
For transmission instruments (suitable for a wide range of various manufacturers) on A105/1 upper basket with A800 central filter.



A813 Adaptor, 19721413
For reprocessing of up to 4 transmission instruments with external spray channel. To be used on A105/1 upper basket.




Adaptor, 06894700, 06895000
For W&H air scaler tips. For W&H Piezo and Prophylaxis tips. Can be connected to A813.



A814 Adaptor, 19721430
For reprocessing of various types of scaler tips. To be used on A105/1 upper basket.



A815 Nozzle, 19721414
For the reprocessing of transmission instruments with an external spray channel. For use with A813 or A105/1 with A800 filter 4 nozzles with silicone hose.

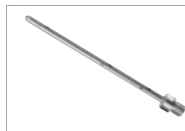

A836 Dental nozzle set, 19721451
Set of 22 injectors for A105/1.



E476 Adaptor, 19721425
Adaptors for use in trays with 5 mm mesh. For transmission instruments with a diameter of 4 to 8 mm (50 per bag).



ADS1 Adaptor silicone, 19721408
Adaptor silicone for AUF1. For connection of transmission instruments with diameter of approx. 20 mm (white).



ADS2 Adaptor silicone, 19721409
Adaptor silicone for AUF1. For connection of transmission instruments with diameter of approx. 16 mm (green).



ADS3 Adaptor silicone, 19721410
Adaptor silicone for AUF1. For connection of transmission instruments with diameter of approx. 22 mm (red).



AUF1 Adaptor, 19721406
Sleeve for transmission instruments (see table at last page).

Filter + Brush



A800 Central filter, 19721404
Filter made from multiply stainless steel mesh for A105/1 upper basket. Mesh size 63 µm. Can be cleaned with A804 cleaning brush. L 405 mm, Ø 22.5 mm



A804 Cleaning brush, 19721412
Twisted brush with spiral bristles for A800 central filter cleaning. Handle with grommet. Length 500 mm, Brush length 100 mm

Inserts



E130 Insert 1/2 for 10 trays, 19721415
11 holders (10 compartments). H 170, spacing 35 mm, H 180, W 180, D 445 mm



E131/1 Insert, 19721416
For optimum loading of 5 mesh and kidney trays. H 168, W 180, D 495 mm



E146 Insert 1/6, 19721417
Mesh size: base: 3 mm, sides: 1,7 mm, lid: 8 mm, 2 hinged handles. For upper or lower basket. H 55, W 150, D 225 mm

Inserts



E197 Insert 1/6 mesh tray, 19721418
Wire mesh. Mesh size: Base 3 mm, sides 1,7 mm, lid 8 mm, H 42, W 150, D 225 mm



E198 insert 1/2, 19721419
For 6 mesh trays/kidney dishes. 7 holders (6 compartments). H 160, spacing 66,5 mm (compatible with mesh tray E197) H 160, W 180, D 495 mm



E337/1 Insert 2/5, 19721420
For instruments in upright position: 4 plastic receptacles, 12 compartments, approx. 22 x 28 mm, 4 compartments, approx. 25 x 28 mm, 48 compartments, approx. 13 x 14 mm H 113, W 173, D 445 mm



E363 Insert 1/6 mesh tray, 19721421
Mesh size 1 mm, with lid (without receptacles for instruments) H 55, W 150, D 225 mm



E379 Insert 1/2 mesh tray, 19721422
Wire mesh with 2 handles. Mesh size 1,7 mm, 80 (+30), W 180, D 445 mm



E430/1 Insert 1/3 mesh tray, 19721423
Wire mesh. Mesh size 5 mm H 40, W 150, D 445 mm



E473/2 Insert, 19721424
Mesh basket with lid for small parts. For hooking onto mesh trays. H 85, W 60, D 60 mm



E520 Insert, 19721426
For 18 root canal instruments hinged. Can be sterilised in steam at 121°C/134°C. H 45, W 75, D 30 mm



E521/2 Insert, 19721427
For 7 extraction or orthopaedic molar forceps. Compartment size 21 x 80 mm. H 135, W 100, D 189 mm



E522/1 Insert, 19721428
9 holders for dental impression plates. H 140, W 100, D 190 mm



E523 Insert 1/2, 19721429
For mesh trays: 7 holders (6 compartments). H 145 Spacing 50 mm, H 151 W 220, D 450 mm



E339/1 Insert, 19721437
Insert for trays

Dispensers



DOS K85/1 200-240V 50Hz, 19721435
Dispenser module



DOS K85/1 NER, 19721442
Dispenser module for Scandinavia



DOS K85/1 200-240V 60Hz, 19721443
Dispenser module

Reprocessing with the ideal equipment package

Upper basket
Flexible for inserts and height adjustable.



1



5

Insert
Intended for different instruments and for flexible use.



2



6

Lower basket
For holding mesh trays as well as various inserts.



3



7-8

Insert
For accommodating various instruments and insertable for upper and lower basket.



4

Insert
Different compartments for optimum loading of upright instruments.

Insert
Equipped with different compartments for optimum loading.

Upper basket – injector nozzles
For various instruments with injector nozzles.

»starter pack« 1-6
19721439
includes the perfect range of basic equipment.

»injection pack« 3-8
19721440
complemented by an additional equipment perfectly and particularly suitable for transmission instruments.

Please note:
The devices are sold empty. Add the appropriate equipment to your device.

Cleaning and disinfection of your transmission instruments

W&H transmission instruments

The special adaptors facilitate an easy cleaning and disinfection and ensure a long lifespan for your transmission instruments.

Following types need additional adaptors:

AM 25 RM	A801 + AUF1 + ADS1 silicone (white)
AM 25 L RM-RM 6-hole	
AM 25 A RM-RM 4-hole	
Alegra TE-97 RM-RM 4-hole	
Alegra TE98 RM-RM 4-hole	
PB-3 (S, LED, LED S), PA-1 LED	A801 + AUF1 + ADS3 silicone (red)
PB-5 (LQ, L, LS), Proxeo Aura PF-60 (Sprayhead 90° and 120° perio short)	A801 + AUF1 + ADS1 silicone (white)



AUF1



ADS1



ADS2



ADS3

Competitors transmission instruments

Clean and disinfect your entire product range efficiently, with or without special adaptors.


Following types need additional adaptors:

RM Connection 4-hole » Turbine, AM	A801 + AUF1 + ADS1 silicone (white)
RM Connection 6-hole » Turbine, AM	
BC 2/3 Hole » Turbine, AM	
NSK Mach®	

No additional adaptor needed:

ISO-Connection (E-Type), Contra-Angles	Only A803 required, no additional adaptor
RQ Connection	
KaVo Multiflex®	
Bien Air Unifix®	
Morita CP4-O	

*NSK Mach®, KaVo Multiflex® and Bien Air Unifix® are trademarks of third parties not affiliated with W&H Dentalwerk Bürmoos GmbH

Please note: Products without the symbol , cannot be reprocessed in the thermal washer disinfectors.

W&H Sterilization Srl.

Via Bolgara, 2
24060 Brusaporto (BG), Italy
t +39 035-66 63 000
f +39 035-50 96 988
office.sterilization@wh.com
wh.com

Pictures are for illustrative purposes only.
© W&H Dentalwerk Bürmoos and © Miele GmbH.
Additional equipment and accessories shown are not included as standard.

20517 AEN Rev. 004 / 19.10.2021
Subject to alterations